

High-Precision, Advanced Tuning Fork Balance

# HJ-E / HJR-E Series

## Operation Manual

### **IMPORTANT**

To ensure safe and proper use of the balance, please read this manual carefully.

After reading this manual, store it in a safe place near the balance, so you can review it as needed.

**VIBRA**

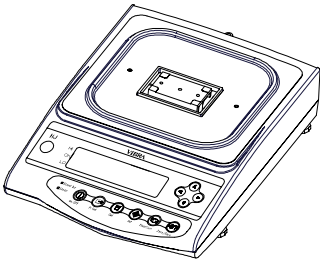
# PREFACE

Thank you for purchasing an HJ-E/HJR-E Series electronic balance.

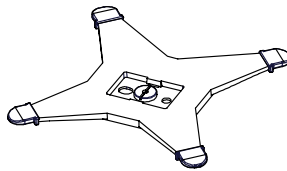
This balance is equipped with a high precision mechanism. It provides parts counting, percentage weighing, and limit functions helpful in weighing a certain amount, and has various other functions. This advanced balance employs a user-friendly program, and the key arrangement is also easy to operate. In addition, the fluorescent display is eye-friendly, and its high-speed, stable performance improves your work efficiency.

## Check for the following accessories before use.

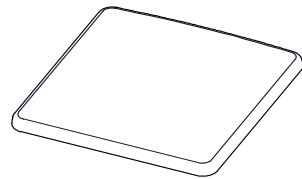
(1) Main Unit



(2) Pan Base

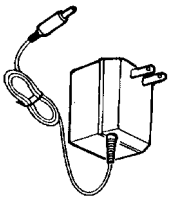


(3) Pan



HJ-21KE/31KE : beforehand installed at factory

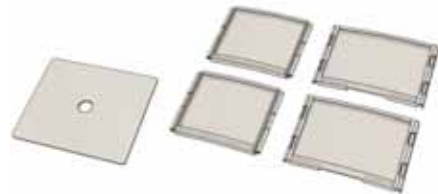
(4) AC adapter



(5) Operation manual



(6) Wind shield panel(5pcs.)



Cover      Panel(2large & 2small)

for 220g ~ 620g Capacity model

# Contents

---

<b>1</b>	<b>Precautions Relating to Use.....</b>	<b>2</b>	<b>9</b>	<b>Calibration and Span Test for the Balance</b>	
			9.1	Span Adjustment.....	37
<b>2</b>	<b>Names of Component Parts .....</b>	<b>5</b>	9.2	Span Test.....	38
			9.3	Calibration of Built-in Weight.....	39
<b>3</b>	<b>Installation of the Balance and Operation Check</b>		<b>10</b>	<b>Date and Time Setup .....</b>	<b>40</b>
3.1	Installation .....	8			
3.2	Windshield setting .....	8	<b>11</b>	<b>Miscellaneous Functions</b>	
3.3	Operation Check .....	9	11.1	Auto Sleep Function.....	41
<b>4</b>	<b>Function 1</b>		11.2	Auto Power Off Function.....	41
4.1	Setting and Check.....	12	11.3	Set Unit Function .....	41
4.2	Description of Function 1.....	13	11.4	Minimum Readability Setting Function..	41
4.3	Additional Functions.....	15	11.5	Advice CAL .....	42
4.4	Interface .....	16	11.6	Date Display .....	42
<b>5</b>	<b>Function 2.....</b>	<b>17</b>	11.7	Time Stamp Output .....	42
			11.8	Save Tare Weight Function.....	42
<b>6</b>	<b>Weighing Mode</b>		11.9	Direct Start Function .....	42
6.1	Measurement Modes.....	18	11.10	Interval Output Function.....	43
6.2	Weighing Machine.....	18	11.11	ID Number Entry.....	44
6.3	Parts Counting .....	19	11.12	Entry of Weight Error .....	45
6.4	Percentage Weighing .....	21	<b>12</b>	<b>Input/Output Functions</b>	
6.5	Unit Converting .....	23	12.1	RS232C Output .....	46
6.6	Gravimeter .....	24	12.2	Output to Peripherals.....	49
<b>7</b>	<b>Addition Function</b>		12.3	Type of Communication Texts.....	50
7.1	Cumulate Function .....	28	12.4	Output Data .....	50
7.2	Net Addition Function .....	29	12.5	Input Commands .....	53
<b>8</b>	<b>Limit Function</b>		<b>13</b>	<b>Use Printers.....</b>	<b>58</b>
8.1	Setting the Limit Function.....	30			
8.2	Judgment and Saving.....	30	<b>14</b>	<b>Output in Compliance with ISO/GLP/GMP .....</b>	<b>59</b>
8.3	Display of Judgment Results.....	30			
8.4	Judge by Absolute Values .....	31	<b>15</b>	<b>Operate on Batteries .....</b>	<b>62</b>
8.5	Judge by Deviation Values .....	33			
8.6	Bar Graph for the 2-point Scale.....	36	<b>16</b>	<b>Troubleshooting.....</b>	<b>63</b>
			<b>17</b>	<b>Specifications .....</b>	<b>64</b>
			<b>18</b>	<b>Conversion Table of Units.....</b>	<b>68</b>

# 1 Precautions Relating to Use

- This Section "Precautions Relating to Use" sets forth precautionary notes that the user should observe in order to prevent physical injury to the user and/or damage to property.
- The nature of problems that may result in the event of improper operation, and consequential effects on the quality and performance of the balance, are indicated under the two categories of "Caution" and "Recommended," and explained using symbols.

## Caution

This is a category to indicate improper handling that may cause physical injury or severe damage to property. Be sure to follow the directions for safe use to avoid serious consequences.

## Recommended

This term indicates steps that the user should take to ensure the quality and reliability of the balance.

### Meanings of Symbols Each symbol is accompanied by an instruction.



Mandatory Symbol

: Indicates a "mandatory" action that should be executed without fail.

Example



Check Level



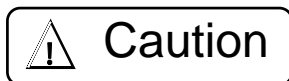
Prohibitive Symbol

: Indicates a "prohibited" action that must not be executed.

Example



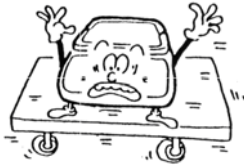
Do Not Use



		<p>◆ <b>Do not disassemble or modify the unit.</b></p> <ul style="list-style-type: none"> <li>• Could cause malfunction or heat generation</li> <li>• Contact our Marketing Division or Technical Service Division.</li> </ul>
		<p>◆ <b>Only AC power should be used.</b>            ◆ <b>Only use the dedicated AC adapter.</b></p> <ul style="list-style-type: none"> <li>• Use of other types of power or adapters may result in heat generation or malfunction of the balance.</li> </ul>
		<p>◆ <b>Do not move the balance when a sample is loaded.</b></p> <ul style="list-style-type: none"> <li>• The loaded sample may fall off the pan and cause an injury.</li> </ul>



Do Not Use



◆ **Do not place the balance on an unstable base or use the balance in a location where it may be subjected to shock.**

- The loaded sample may fall off the pan.
- Accurate measurement may be rendered impossible.



Do Not Drop



◆ **Do not lay the AC adapter cable on the surface of the passage.**

- Somebody may trip on the cable, causing the balance to fall off, thereby causing injury and/or damage to the balance.



Do Not Leave Afloat



◆ **Do not use the balance with its adjusters lifted.**

- The balance will become unstable, preventing accurate measurement.



Don't touch with wet hand



**Do not touch the AC adapter or balance with wet hands.**

- Danger of electric shock.



Don't make it wet



**Do not use the balance in a location where it may be subjected to excess moisture.**

- Electric shock or short-circuiting could occur.
- The balance may be corroded, with resultant malfunction.



Avoid excess dust



**Do not use the balance in a location where it may be subjected to excess dust.**

- Risk of explosion or fire.
- Short-circuit or lack of continuity may occur, leading to a malfunction of the balance.

## Recommended



Calibrate Balance



◆ **Calibrate the balance after installation or relocation.**

- Measurement values may contain errors, preventing accurate measurement from being conducted.



Do Not Apply Force



◆ **Avoid applying excess force or impact to the balance.**

- Place the sample to be measured on the balance carefully to prevent breakage or malfunction.



Do Not Use



◆ **Do not use the balance in a location where it may be subjected to abrupt changes in ambient temperature or humidity.**

- Accurate measurement may be rendered impossible.
- Use the balance in an ambient temperature range of 5°C to 35°C and with 80% or lower relative humidity.



Do Not Overload



◆ **Do not leave the balance overloaded (When it is overloaded, (E r r) is displayed).**

- Take down the loaded sample immediately to prevent breakage or malfunction.



Do Not Use



◆ **Do not use the balance in a location where it is subject to direct sunlight.**

- The indications would be illegible.
- An internal temperature increase in the balance may lead to inaccurate measurement.



Unplug Adapter



◆ **If the balance is to be unused for an extended period of time, unplug the adapter.**

- This conserves power and prevents deterioration.



Do Not Use

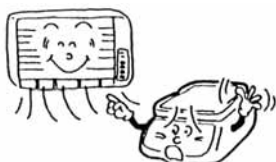


◆ **Do not use volatile solvents for cleaning.**

- The body may be distorted.
- To clean the unit of stains, use a piece of dry cloth or cloth soaked in a small quantity of neutral detergent.



Do Not Use



◆ **Do not use the balance in a location where it may be subject to air from an air-conditioning unit.**

- Extreme changes in the ambient temperature may result in inaccurate measurements.



Do Not Use



◆ **Do not use the balance on a soft floor.**

- When loaded with a sample, the balance may tip or move, preventing accurate measurements from being conducted.



Check Level



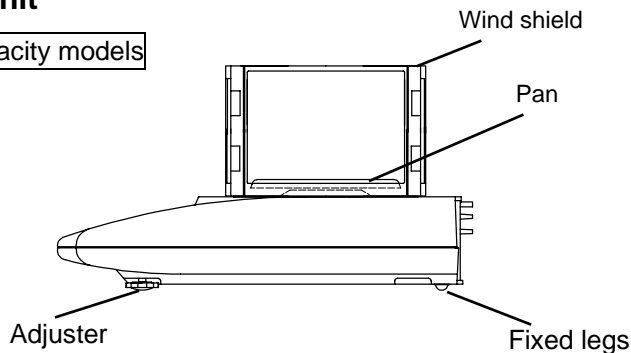
◆ **Do not use the balance when it is tilted.**

- When the balance is tilted, an error may be caused, preventing accurate measurement from being conducted. Place the balance on a level surface.

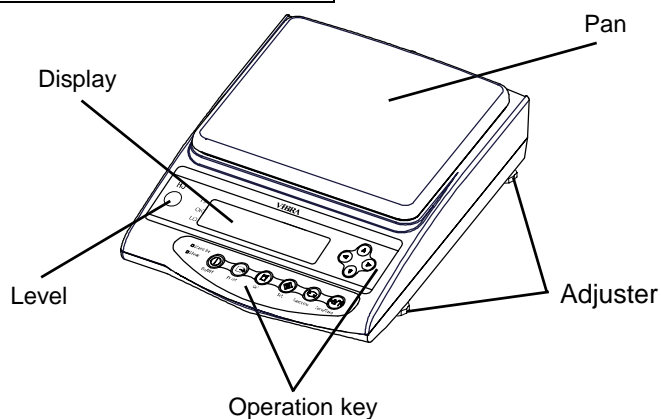
## 2 Names of Component Parts

### 2.1 Main Unit

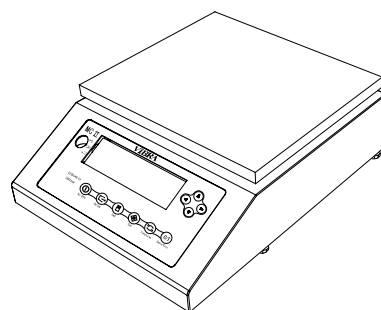
220g ~ 620g capacity models



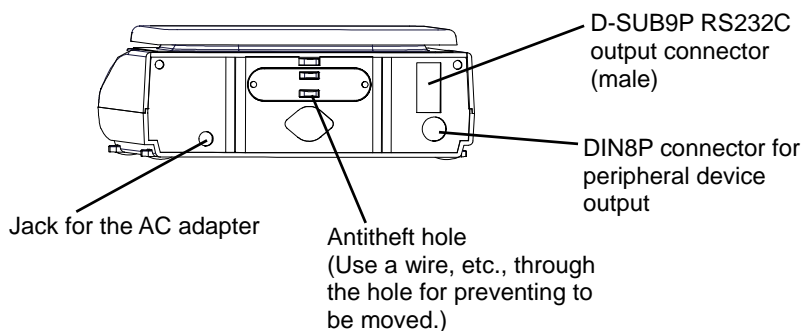
1200g ~ 15kg capacity models



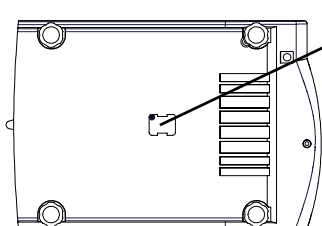
21kg,31kg capacity model



Rear



Bottom



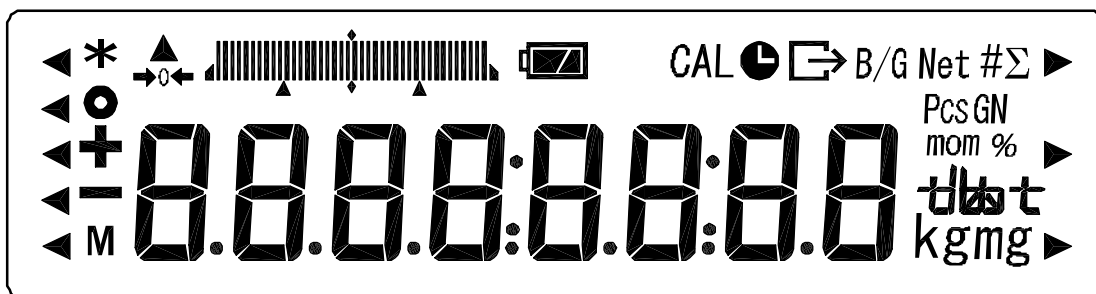
Cover for a hanging hole.

Inside is installed hook for under weighing for 220g ~ 620g and 21kg,31kg model.

For 1200g ~ 15kg capacity model, optional under weighing hook is necessary

\* It is necessary to perform Span adjustment in the case of measuring with a lower plate using a lower plate. When NET is displayed, weighing range is reduced and we cannot weigh up to the full capacity. In this case, use external weight having 50% or more of the full capacity and less enough not to show " E r r " for such adjustment or test.

## 2.2 Displayed Signs and Operation Keys

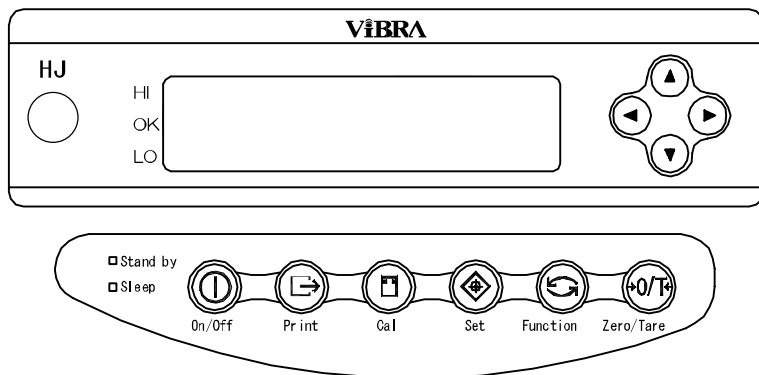









### 2.2.1 Displayed signs

Display	Description
kg g	Kilogram, gram
→0←	Zero-point
-	Minus
Net	Tare sign
B/G	Gross weight
○	Indication of stable balance (If this light is off, the balance is unstable.)
*	Indicates that the addition function is enabled for accepting an additional load.
Pcs	Parts counting
GN	Grain unit
%	Percentage weighing
#	Unit converting
Σ	Displays a sum total.
↔	Displayed when data is output in compliance with ISO/GLP/GMP.
🕒	Lights up when date/time is being set or displayed. Blinks during interval output.
◀	Judgment result display when the limit function is used
mom	Momme
M	Indicates that the value is set in memory (If it is flashing, the value is being saved.)
CAL	Lights or blinks upon span adjustment or when Advice CAL is in operation.
	<ul style="list-style-type: none"> <li>Bar graph (Refer to Hints on page 36).</li> </ul>
	Displays a set unit.
	Battery indicator. As the battery gets weaker, the sign is changed into [  ] or [  ]. (Refer to Section 15: Operate on Batteries on page 62).



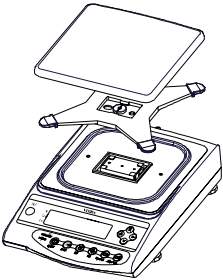
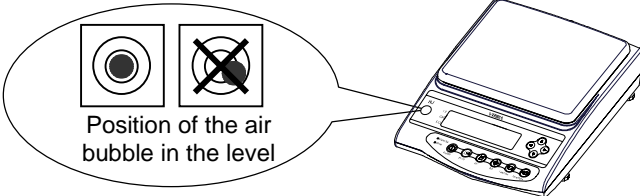
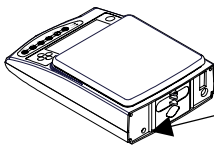
## 2.2.2 Operation keys and their functions



Operation key		Function
	<b>On/Off</b> key	Turns the balance on or off.
	<b>Print</b> key	Starts output. Also used to cancel the operation to set date/time.
	<b>Set</b> key	[Short press] Starts the setting a count or a percent. [Continuous press] Starts the setting of a limit value when the limit function is enabled.
	<b>Function</b> key	[Short press] Switches the measurement mode. [Short press] Used to enter numerical values. [Short press] Used to select a function to set. [Continuous press] Calls functions.
	<b>Zero/Tare</b> key	[Short press] Used for zero-setting or sets the display to zero by tare range. [Short press] Used to enter numerical values. [Short press] Used to select a function.
	<b>Cal</b> key	Starts span adjustment or a span test.
	Arrow keys	The arrow keys function in the same way as the <b>Function</b> key or the <b>Zero/Tare</b> key when you set a function or enter numerical values for additional functions.
LED (green)	Stand by	Lights up when the power is off (standby). * The LED does not light up when the balance is operated on batteries.
LED (orange)	Sleep	Lights up when the auto sleep function is enabled.

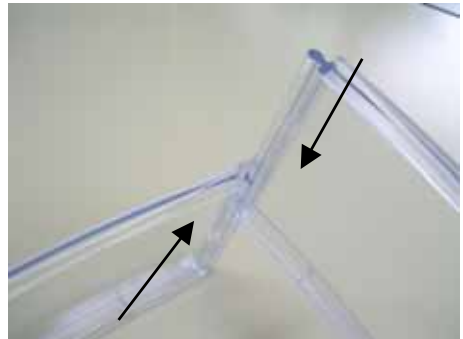
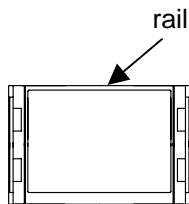
# 3 Installation of the Balance and Operation Check

## 3.1 Installation

<p>1 Attach the pan base and pan to the main unit.</p> 	<p>Attach the pan base to the main unit. Be sure that the "FRONT" indication is on the display side when attaching the pan base. Align the bosses on the body with the holes on the pan base before fastening the screws. After attaching the pan base, put the pan on it.</p>
<p>2 Level the balance.</p>  <p>Position of the air bubble in the level</p>	<p>Turn the adjusters so that the air bubble of the level is within the red circle.</p> <p>&lt;Caution&gt; Be sure that none of the adjusters is lifted from the floor.</p>
<p>3 Connect the AC adapter.</p>  <p>Jack for the AC adapter</p>	<p>Connect the AC adapter to the balance.</p> <p>&lt;Caution&gt; (Refer to Section 15, "Operate on Batteries," on page 62 if the battery option is employed.)</p>

## 3.2 Windshield setting ( for 220g ~ 620g capacity models )

1. On to the rail, put the panels




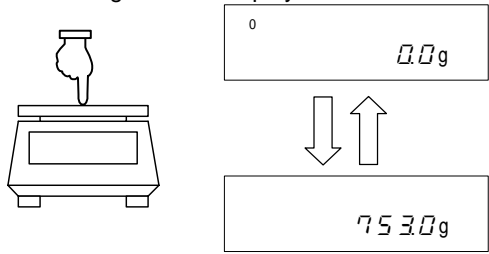

2. After putting 4 panels in the same way, put the cover

3. Place the windshield along the guide.

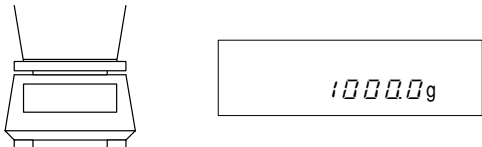
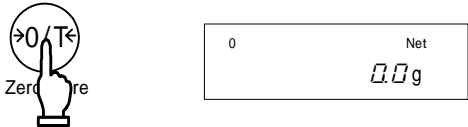
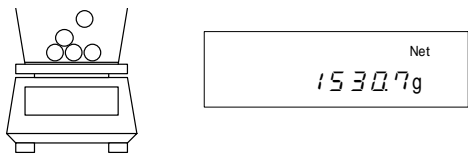


## 3.2 Operation Check

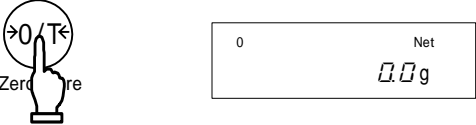
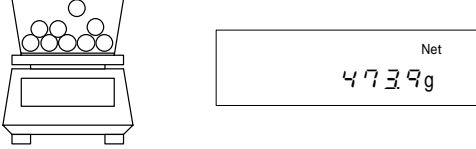
### 3.2.1 Turn the power on and off and check the display

	<p>Connect the AC adapter. The balance enters standby mode, and the Stand by lamp (LED) lights up.</p> <p>Press the <b>On/Off</b> key. When the display is turned on, check it for any missing signs.</p>
<p>Check changes of the display.</p> 	<p>Give the pan a slight press, and check that the display is changed.</p> <p>Check that the display goes back to zero when you lift your hand.</p>
<p>Stand by</p> 	<p>Press the <b>On/Off</b> key again. The balance enters standby mode, and the Stand by lamp (LED) lights up.</p>

### 3.2.2 Set a Tare Weight

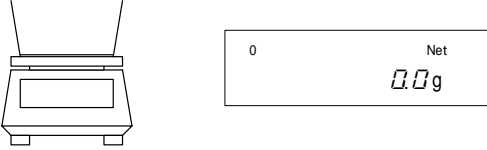
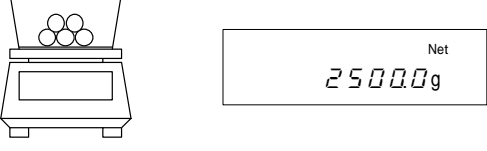
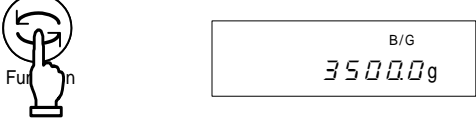
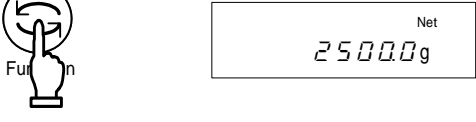
<p>1 Put tare (container) on the balance.</p> 	<p>Put tare on the balance, and the weight is displayed.</p>
<p>2 Reset the display to 0.</p> 	<p>Press the <b>Zero/Tare</b> key. The tare weight is set, and the display is reset to zero.</p>
<p>3 Put a sample on the balance.</p> 	<p>The net weight of the sample in the container is displayed.</p>

## ☆ Weigh additional samples

<p>4 Reset the display to 0.</p> 	<p>Press the <b>Zero/Tare</b> key. The gross weight of the sample on the pan is reset to 0.</p>
<p>5 Put an additional sample on the balance.</p> 	<p>The weight of the additional sample is displayed.  In this way, you can add a sample and weigh it by pressing the <b>Zero/Tare</b> key to reset the displayed weight to 0.</p>

### 3.2.3 Display a gross weight

This is the function to display gross weight of samples with tare. Sample weight with tare is referred to as gross weight, and sample weight excluding tare is referred to as net weight. Gross weight can only be displayed when the balance is used as a weighing machine.

	<p>Put tare on the balance to set the tare weight.</p>
	<p>Put a sample to weigh on the balance. At this point, the net weight of the sample is displayed.</p>
	<p>Press the <b>Function</b> key once.  The gross weight including the tare weight is displayed. The B/G sign is displayed.</p>
	<p>When you press the <b>Function</b> key again, the balance is reset to the normal mode.</p>

### Caution:

1. Gross weight can only be displayed when the balance functions as a weighing machine.
2. While gross weight is displayed, you cannot set the tare. You can only adjust the zero-point.

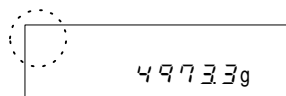
☆ Hints ☆

The following description applies commonly to the weighing machine, parts counting, percentage weighing, unit converting, and gravimeter functions.

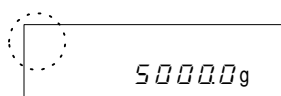
1. A small amount of current is still flowing after the balance is turned off.  
In a standby state, the Stand by lamp (LED) is lit. When the balance is turned on, the LED goes out.  
\* The LED does not light up when the balance is operated on batteries.
2. The bar graph indicates the current load condition in relation to the weighing capacity. As the load approaches the weighing capacity, the bar gets closer to the right end.  
\* When a tare range is set, the weight is still reflected in the bar graph even if the display indicates 0.
3. When the balance is stable, a circle is displayed. When it gets unstable, the sign goes out. When the balance is affected by wind or vibration, numeric values on the display may flicker or the "stable" sign may go out. In such a case, reset Function 1 (refer to page 13) to improve the stability by referring to the following table.

Influence of wind or vibration	Items of Function 1		
	Stability Judgment 4 5.d	Response Speed 5 r E	Setting the lowest digit of display for Unit A (B) [b 2, dR] [b 4, db]
Small	1	1	1
↕	2	2	2
↕	3	3	3
↕	4	4	4
Large		5	5

If the influence of the wind or vibration is small, set these setting items of Function 1 higher. Set these items lower as the influence is increased.



Unstable



Stable

4. When you reset the display to 0 or set a tare range, the display indicates 0, and the [→0←] sign is displayed. When you set a tare range, the [Net] sign is displayed.



\* When the zero-point deviates from the true zero-point by one fourth of a division, the →0← sign goes out






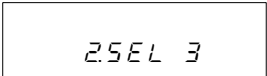



\* When you set a tare range, the display indicates 0, and the Net sign is displayed.

5. When a tare range is set, the weighing capacity is reduced accordingly. Weighing capacity = original weighing capacity – tare weight
6. When [d - E r r] is displayed after an object is put on the balance, it indicates that the weighing capacity is exceeded.
7. When you remove the pan, and the zero-point falls below the original zero-point, [u - E r r] is displayed.
8. With any of the parts counting function, percentage weighing function, or unit converting function, the display still indicates zero when you press the pan unless any weight data of a sample is stored in the memory.
9. When you turn on the balance, it starts in the mode when it was turned off. For example, if you turn off the balance in the parts counting mode, it starts in the parts counting mode when it is turned on the next time.

# 4 Function 1

## 4.1 Setting and Check

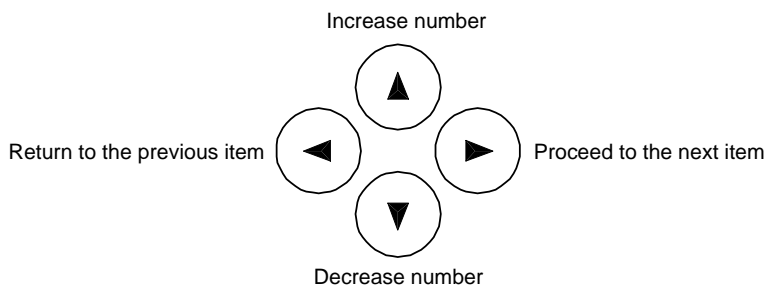
<p>1. Call Function 1.</p>  <p>Continuous press</p>  <p>Key released</p>	<p>Keep pressing the <b>Function</b> key. When the display is changed to [ <b>Func</b> ], release the key.</p> <p>Now you can set Function 1. The first setting item is displayed [ <b>1 SET 1</b> ]. (Refer to Section 4.2: Description of Function 1 on page 13.)</p>
<p>2. Select the next setting item.</p>  	<p>Every time you press the <b>Function</b> key, you are moved one item forward.</p>
<p>3. Change settings.</p>  	<p>Every time you press the <b>Zero/Tare</b> key, the value in the right end is changed. Make appropriate settings.</p>
<p>4. Complete the setting of functions.</p> 	<p>Press the <b>Set</b> key, or press the <b>Function</b> key several times until the balance enters the measurement mode.</p> <p>The setting of functions is completed, and the balance goes back into the measurement mode.</p>

\* The setting values of Function 1 are held after the power is turned off.

### Function setting using the arrow keys

Using the arrow keys located on the right side of the display can make function setting operations easier.

You can use the arrow keys instead of the **Zero/Tare** key or the **Function** key when any Function 1 item is displayed. The arrow keys can be used as shown in the following figure:



## 4.2 Description of Function 1

Item		Set Value	Description		
Weighing Mode		1. <i>SEt.</i>	☆1	Weighing machine (only supports weight measuring)	
			2	Parts counting (parts counting and weight measuring)	
			3	Percentage weighing (weight percent measuring and weight measuring)	
			4	Unit converting (coefficient multiplying and weight measuring)	
			5	Gravimeter (measurement of specific gravity)	
Displayed for the gravimeter.	Media	11. <i>NEd.</i>	☆0	Water	
			1	Any liquid other than water	
	Output data	12. <i>d.o.d.</i>	☆0	Only a specific gravity value	
			1	Specific gravity, weight, and water temperature or the specific gravity of the media (liquid)	
	Auto output	13. <i>RO.</i>	☆0	Disabled (Print key is used for output)	
			1	Output once after a specific gravity measurement.	
Additional Functions		2. <i>SEL</i>	☆0	Disable additional functions.	
			1	Addition function [ <i>2C. RdR.</i> ]	
			2	Limit function	
			3	Addition and limit functions	
Auto-Zero (Zero Tracking)		3. <i>RO</i>	0	Disable	Functions to automatically adjust for a slight zero-point deviation.
			☆1	Enable	
Stability Judgment		4. <i>S.d.</i>	1	Wide (Mild)	
			☆2	↓	
			3	Narrow (Strict)	
			4		
Response Speed		5. <i>rE.</i>	0	Measurement by consecutive weighings	
			1		
			2	Fast	
			☆3	↓	
			4	Slow	
Interface		6. <i>IF.</i>	0	Stop input/output	
			☆1	Numeric 6-digit format	
			2	Numeric 7-digit format	
			3	Extended 7-digit numeric format	
Span Adjustment Span Test		7. <i>CR.</i>	0	Disable the Cal key.	
			☆1	1	Span adjustment with built-in weight (Auto-span adjustment)
			2	Span test with built-in weight	
			☆2	3	Span adjustment with external weight
			4	Span test with external weight	Displayed only on the HJR series balance.
Bar Graph		8. <i>b.G.</i>	0	No display	
			☆1	Displays the bar graph	
Auto Power Off		9. <i>RP.</i>	0	Disable (for continuous use)	Available only when operated on batteries.
			☆1	1	
Auto sleep		10. <i>RS.</i>	0	Disable	
			☆1	1	The balance goes into sleep mode 3 minutes later when it is operated on AC power.

A star (☆) denotes a factory setting.

☆1 denotes a factory setting for the HJR series, and ☆2 denotes a factory setting for the HJ series.

Item	Set Value	Description	
Unit A	<i>b 1 uR</i>	☆ 1	[g]
		2	[kg]
		4	[ct] (ct)
		5	[oz] (oz)
		6	[lb] (lb)
		7	[oz t] (ozt)
		8	[dwt] (dwt)
		9	[GN] (grain)
		R	[tl] (tl_Hong Kong)
		b	[tl ▶ Upper right] (tl_Singapore, Malaysia)
		L	[tl ▶ Lower right] (tl_Taiwan)
		d	[mom]
		E	[to] (to)
		Unit A Min. indication setting	<i>b 2 dR</i>
2	↓		
3	Rough		
4			
5			
Unit B	<i>b 3 u.b</i>	☆ 0	None
		1	[g]
		2	[kg]
		4	[ct] (ct)
		5	[oz] (oz)
		6	[lb] (lb)
		7	[oz t] (ozt)
		8	[dwt] (dwt)
		9	[GN] (grain)
		R	[tl] (tl_Hong Kong)
		b	[tl ▶ Upper right] (tl_Singapore, Malaysia)
		L	[tl ▶ Lower right] (tl_Taiwan)
		d	[mom]
		E	[to] (to)
Unit B Min. indication setting	<i>b 4 d.b</i>	☆ 1	Fine
		2	↓
		3	Rough
		4	
		5	
		* To be displayed when [ <i>b 3 u.b</i> ] is other than [0].	
Advice CAL	<i>d Rd.C</i>	☆1 0	Disable
		☆2 1	Enable
Compliance with ISO/GLP/GMP	<i>E GLP</i>	☆ 0	Disable
		1	Enable

A star (☆) denotes a factory setting.

☆1 denotes a factory setting for the HJ series, and ☆2 denotes a factory setting for the HJR series.



Displayed when [E.L.P.] is [1].	Output of CAL and span test results	E1 out	0	Disable
			★1	Enable
	Data compliant with GLP	E2 od	★0	Disable
			1	Enable
Printed Language	E3 P.F.	★1	English	
		2	Japanese (Katakana)	
Date Display	F. dRtE	1	Output in Year-Month-Day format.	
		2	Output in Month-Day-Year format.	
		★3	Output in Day-Month-Year format.	
Time Stamp Output	G. t.o.	★0	Disable	
		1	Outputs time together with measurement data.	
Set Tare Range	H. tR	1	Reset the display to 0 immediately after the Zero/Tare key is pressed.	
		★2	Reset the display to 0 when the balance is stabilized after the Zero/Tare key is pressed.	
Store Tare	J. tRrE	★0	Disable	
		1	Enable	
Direct Start	L. dSt.	★0	The balance goes into the standby mode when the AC adapter is plugged in.	
		1	The balance is turned on when the AC adapter is plugged in.	

### 4.3 Additional Functions

Displayed when [2. SEL] is other than [0].

Item	Set Value		Description
Condition	21 Ca	★1	Always judge (even when the balance is unstable).
		2	Judge only when the balance is stable.
Range to Cover	22 L i	0	Above +5 ticks is covered in the judgment. (The range +5 ticks or less is not covered for judgment including negative values.)
		★1	Full range judgment (including negative values)
Point Scale	23 P i	1	1-point scale (OK/LO)
		★2	2-point scale (HI/OK/LO)
		3	3-point scale (scale of 1 to 4)
		4	4-point scale (scale of 1 to 5)
Judge by	24 tYP.	★1	Judge by absolute values.
		2	Judge by deviation values.
Buzzer for rank 1	25 bu.1	★0	Buzzer is turned off for rank 1 (LO).
		1	Buzzer turns on for rank 1 (LO).
Buzzer for rank 2	26 bu.2	★0	Buzzer is turned off for rank 2 (OK).
		1	Buzzer turns on for rank 2 (OK).
Buzzer for rank 3	27 bu.3	★0	Buzzer is turned off for rank 3 (HI).
		1	Buzzer turns on for rank 3 (HI).
Buzzer for rank 4	28 bu.4	★0	Buzzer is turned off for rank 4.
		1	Buzzer turns on for rank 4.
Buzzer for rank 5	29 bu.5	★0	Buzzer is turned off for rank 5.
		1	Buzzer turns on for rank 5.

How to Indicate Results	2R LG	☆ 1	Displays HI, OK, LO, or the rank.
		2	Bar graph for the 2-point scale (Available only when you select the 2-point scale)
Relay output control <sup>*1</sup>	2b r.o.c.	☆ 1	Output all the time (independently of external signal input).
		2	Controlled by an external input signal.
Addition operation <sup>*2</sup>	2C AdA	☆ 1	Cumulate function
		2	Net addition function

\*1 Set [2b r.o.c. 2] when the limit contact output option or the full pack option is used. Otherwise, use [2b r.o.c. 1].

\*2 Displayed when [2 SEL] is [1] or [3].

## 4.4 Interface

Displayed when [5 IF] is [1~3].

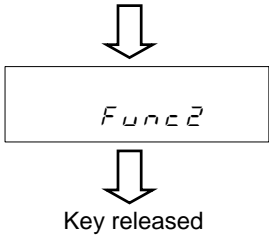
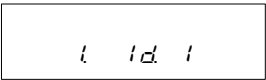
Item	Set Value	Description	
Output Control	5I o.c.	0	Stop output
		1	Output continuously at all times
		2	Output continuously if stable (Stop output if unstable)
		3	Output once when the <b>[Print]</b> key is pressed (regardless of whether the balance is stable or unstable).
		4	Output once when the balance is stable. After any sample is lifted causing the display to show zero or less value, outputs the next time any sample is put on again and the balance is stabilized.
		5	Output once when the balance is stable. Stop output when unstable. Outputs once when the balance is stabilized again (including the zero point) even if it is not reloaded.
		6	Output once when the balance is stable. Output continuously when unstable. Output is stopped after a single output when the balance is stable even if it is not reloaded.
		☆7	Output once when the <b>[Print]</b> key is pressed if the balance is stable.
		A	Output once immediately after a specified interval.
		b	Output once immediately after a specified interval if the balance is stable.
Baud Rate	52 b.L.	☆ 1	1200 bps
		2	2400 bps
		3	4800 bps
		4	9600 bps
		5	19200 bps
Parity	53 P.R.	☆ 0	None
		1	Odd
		2	Even
Data Length	54 d.L.	7	7 bits
		☆ 8	8 bits
Stop Bits	55 St.	1	1 bit
		☆ 2	2 bits
Unused Higher Digits	56 u.r.	☆ 0	Embed 0 (30H).
		1	Filled by space characters (20H).
Response Commands	57 r.E.S.	☆ 1	The A00 and Exx responses are used.
		2	The ACK and NAK responses are used.

A star (☆) denotes a factory setting.

\* The data output interval is 0.1 to 1 second for continuous output. (The interval varies with the measurement conditions, etc.)

## 5 Function 2

### 5.1 Setting and Check

<p>1 Call Function 2.</p> <p>Press the <b>Function</b> key while holding down the <b>Zero/Tare</b> key.</p> 	<p>Press the <b>Function</b> key while holding down the <b>Zero/Tare</b> key.</p> <p>When [<i>Func 2</i>] is displayed, release the key.</p> <p>The display will be changed to show [<i>t id 0</i>].</p>
<p>2 Change settings.</p> 	<p>Use the <b>Function</b> key to select a setting item.</p> <p>Change the setting in the right end with the <b>Zero/Tare</b> key.</p>
<p>3 Complete the setting of functions.</p>	<p>Press the <b>Set</b> key, or press the <b>Function</b> key several times until the balance enters measurement mode.</p> <p>The setting of Function 2 is completed, and the balance goes back into the measurement mode.</p>

### 5.2 Description of Function 2

Item	Set Value	Description		
* Setup of ID No.	<i>t id</i>	☆0	Disable	
		1	Enable	
* Setting instrumental error of the weight in use	<i>2. o.n.P.</i>	☆0	Disable	
		1	Enable	
* Calibration of Built-in Weight	<i>3. r.C.R.</i>	☆0	Disable	Only supported by the HJR series
		1	Enable	
Use of weight error	<i>4. n.E.H.</i>	☆0	Do not use.	
		1	Weight error values set are used upon span / span adjustment to use external weights or when a built-in weight is to be calibrated.	

A star (☆) denotes a factory setting.

#### **Caution:**

The setting values for the items marked by an asterisk (\*) are reset to factory defaults every time you call it.

## 6 Weighing Mode

You can select a weighing mode from the five types in Function 1: a weighing machine, a parts counting, a percentage weighing, a unit converting, and a gravimeter. All the weighing modes except for the gravimeter support additional functions (the addition function and the limit function) depending on the display unit.

### 6.1 Measurement Modes

When you select a weighing mode and additional functions, you can weigh samples in different ways as shown below. Every time you press the **Function** key, the measurement mode is switched.

Weighing Mode	Measurement Mode		Additional Functions		Remarks
	Unit	Function	Addition	Limit	
Weighing machine	Unit A	Weight measuring	○	○	
	Unit A (B/G)	Gross weight	×	×	
	Unit B	Weight measuring	×	×	
	Unit A ( $\Sigma$ )	Cumulative weight	Display	×	
Parts counting	Pcs	Counting	○	○	
	$\Sigma$ Pcs	Cumulative count	Display	×	
	Unit A (Pcs)	Average unit weight	×	×	
	Unit A	Weight measuring	×	×	
Percentage Balance	%	Percent measuring	○	○	
	$\Sigma$ %	Cumulative percent	Display	×	
	Unit A	Weight measuring	×	×	
Unit converting	#	Coefficient multiplying	○	○	
	$\Sigma$ #	Cumulative sum	Display	×	
	Unit A	Weight measuring	×	×	
Gravimeter	g	Measurement of specific gravity	×	×	The weight unit is fixed to "gram."

\* ○: Supports both measurement and display.

×: Does not support measurement or display.

Display: The addition value is displayed.

You can select g, kg, or another unit for Unit A and Unit B.

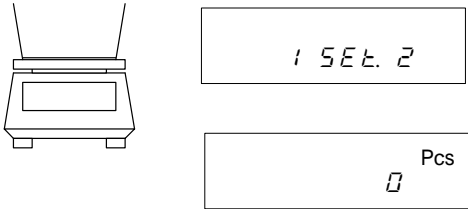
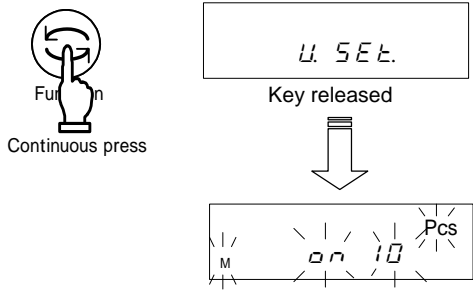
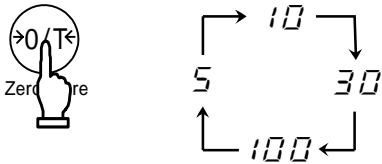
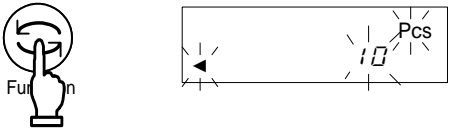
### 6.2 Weighing Machine

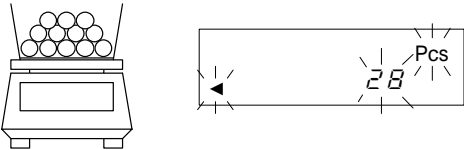
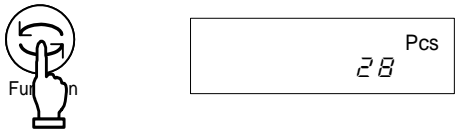
The weighing machine function only supports g, kg, and other units of weight. This function also allows the display of a gross weight (refer to page 10).

## 6.3 Parts Counting

The parts counting function employs an automatic memory update method, also referred to as a Self Counting System. When you put a specified number of samples on the balance, and put additional samples, the balance automatically updates the average unit weight of the samples. However, you cannot add more than three times as many samples as the first reading. This mechanism allows accurate counting.

### 6.3.1 Sampling

<p>1 Select the parts counting function.</p> 	<p>Select the parts counting function [ 1 SEt. 2 ] in Function 1.</p> <p>Then, press the <b>Set</b> key to return to the measurement mode.</p> <p>The weighing mode is changed to "parts counting". Place a tare container on the balance as appropriate.</p>
<p>2 Start sampling.</p> 	<p>Press the <b>Function</b> key for a few seconds.</p> <p>When the display indicates [ u. SEt. ], release the key.</p> <p>The display then blinks [ 0n 10 Pcs ]. This indicates to put 10 samples on the balance.</p>
<p>3 Change the number of samples as required.</p> <p>How the value changes</p> 	<p>If weight variation is significant between samples or unit weight is light, press the <b>Zero/Tare</b> key to sufficiently increase the number of samples.</p> <p>Each time you press the key, the value in the right end is changed. Select an appropriate value.</p> <p>If you do not need to change the number, skip this step.</p>
<p>4 Put a sample on the balance.</p> 	<p>Put the displayed number of samples on the balance, and press the <b>Function</b> key.</p> <p>The display then blinks [ 10 Pcs ].</p>

<p>5 Put additional samples on the balance.</p> 	<p>Add no more than three times as many samples as the displayed number and wait until it is stable again. A beep will sound and the memory is updated. Repeating this operation will improve the resolution of parts counting, resulting in more accurate measurements.</p>
<p>5 Save the average unit weight of the samples.</p> 	<p>Press the <b>Function</b> key to complete the update of the memory. The balance beeps to indicate that the average unit weight is saved, and then goes back into the measurement mode.</p>

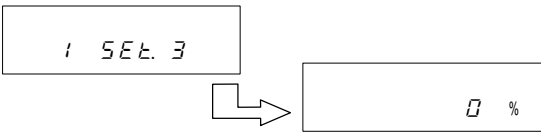
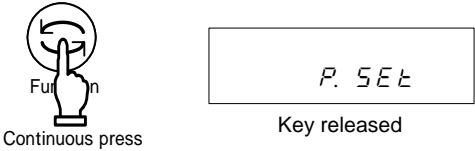
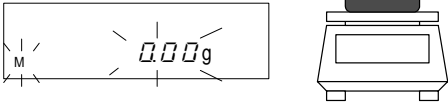

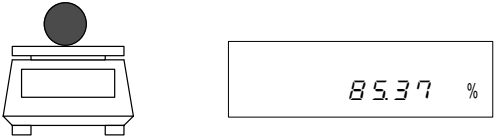
### ☆ Hints ☆

1. The [*S u b*] display that may be shown during a memory update of average unit weight indicates that you added more than three times as many samples as the displayed number and thus the counting accuracy is poor. Do not add too many samples at a time. Gradually increase the number of samples to increase counting accuracy.
2. The [*R d d*] display during a memory update of average unit weight indicates that the number of samples you added is too small and the counting accuracy is poor. Add samples until the [*R d d*] display is turned off to improve the counting accuracy.
3. Even when these indications are displayed, you can press the Function key at any time to store the average unit weight. In this case, however, counting accuracy is low, and a counting error may be caused.
4. When [*L - E r r*] is displayed, it indicates that the average unit weight is too small for the balance to count the samples. For more information on the minimum unit weight supported by the parts counting function, refer to Section 17: Specifications on page 64.

## 6.4 Percentage Weighing

### 6.4.1 Set a reference weight by weighing an actual sample

You can make the display indicate the weight percentage (%) of a sample with respect to a reference sample weight.

<p>1 Select the percentage weighing function.</p> 	<p>Select the percentage weighing function [1 SEt 3] in Function 1. Then, press the <b>Set</b> key to return to the measurement mode. The weighing mode is changed to "percentage weighing". Place a tare container on the balance as appropriate.</p>
<p>2 Set a reference weight.</p> 	<p>Press the <b>Function</b> key for a few seconds. When the display indicates [P. SEt], release the key.</p>
<p>3 Put a sample on the balance.</p> 	<p>The weight display begins to blink. If you have already set a reference weight, that value is displayed. Put a reference sample on the balance.</p>
<p>4 Save the reference weight.</p> 	<p>Press the <b>Function</b> key. The balance beeps to indicate that the value was saved, and then goes back into the measurement mode.</p>
<p>5 Put a sample to weigh on the balance.</p> 	<p>The display indicates the weight percentage (%) of the sample with respect to the reference sample.</p>

### Caution:

The least significant digit of the percent value is automatically selected in accordance with the reference weight saved.

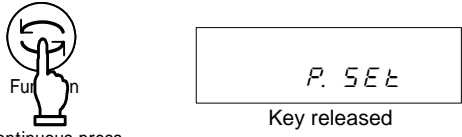
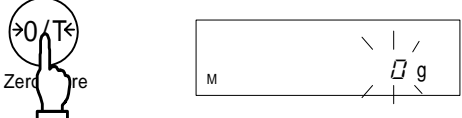
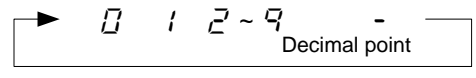
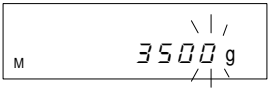

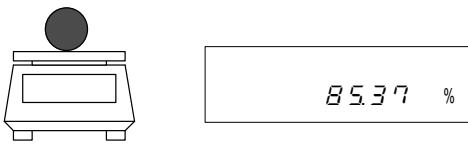
Min. display digit	Range of Reference Weight
L - E r r	Reference weight < Lower weight limit
1%	Lower weight limit ≤ Reference weight < Lower weight limit × 10
0.1%	Lower weight limit × 10 ≤ Reference weight < Lower weight limit × 100
0.01%	Lower weight limit × 100 ≤ Reference weight

You cannot set a value below the lower weight limit as a reference weight.

\* The lower weight limit depends on the model. For more information, refer to Section 17: Specifications on page 64.

## 6.4.2 Set a reference weight by entering a value

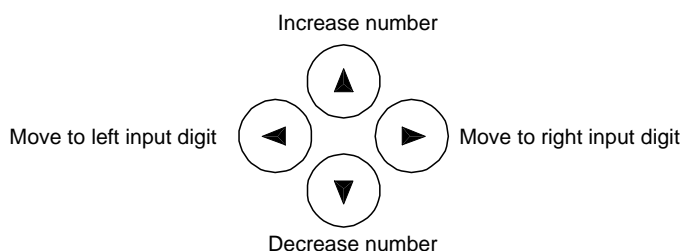
You can enter a reference weight in the balance, and then make the display indicate the weight percentage (%) of a sample with respect to the reference.

<p>1 Set a reference weight.</p>  <p>Continuous press</p>	<p>Press the <b>Function</b> key for a few seconds. When the display indicates [P. S E t.], release the key.</p> <p>The weight display begins to blink. If you have already set a reference weight, that value is displayed.</p>
<p>2 Enter a reference weight.</p> 	<p>Press the <b>Zero/Tare</b> key. A (0) blinks in the right end of the display.</p>
<p>3 Select a number.</p> 	<p>Every time you press the <b>Zero/Tare</b> key, the number changes.</p>
<p>4 Digit selections to enter a reference weight.</p> 	<p>Press the <b>Function</b> key so the number entered is shifted to the left, and you can enter the next digit.</p>
<p>5 Save the reference weight.</p> 	<p>Press the <b>Set</b> key. The display indicates the weight percentage (%) of the sample with respect to the reference weight.</p>
<p>6 Put a sample to weigh on the balance.</p> 	<p>Put a sample to weigh on the balance. The display indicates the weight percentage (%) of the sample with respect to the reference.</p>

### Use the arrow keys to enter values

You can use the arrow keys to enter values instead of the **Zero/Tare** key or the **Function** key.

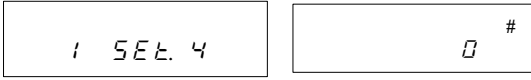
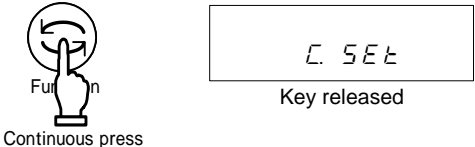

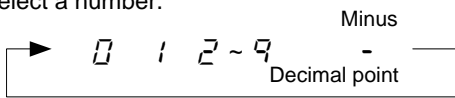

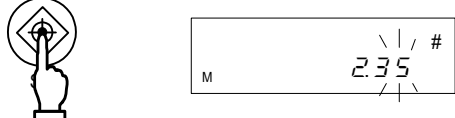
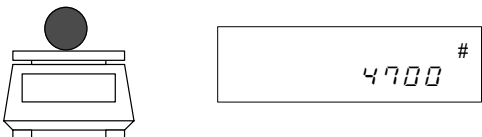
The arrow keys can be used as shown in the following figure:





## 6.5 Unit Converting

The unit converting function multiplies the weight of a sample on the balance by a saved coefficient, and displays the result.

<p>1 Select the unit converting function.</p> 	<p>Select the unit converting function in Function 1 [1 SEt. 4]. Then, press the <b>Set</b> key to return to the measurement mode.</p> <p>The weighing mode is changed to "unit converting." Place a tare container on the balance as appropriate.</p>
<p>2 Set a coefficient.</p> 	<p>Press the <b>Function</b> key for a few seconds.</p> <p>When the display indicates [C. SEt.], release the key.</p> <p>The display blinks a coefficient. If you already have set any coefficient, that value is displayed.</p>
<p>3 Enter a weight.</p> 	<p>Press the <b>Zero/Tare</b> key.</p> <p>A [0] blinks in the right end of the display.</p>
<p>4 Select a number.</p> 	<p>Every time you press the <b>Zero/Tare</b> key, the number changes.</p>
<p>5 Select digit places to specify a coefficient.</p> 	<p>Press the <b>Function</b> key so the number entered shifts to the left and you can enter the next digit.</p>
<p>6 Save the coefficient.</p> 	<p>Press the <b>Set</b> key.</p> <p>The balance beeps to indicate that the value was saved, and then goes back into the measurement mode.</p>
<p>7 Put a sample to weigh on the balance.</p>  <p>Example) 4700 (displayed) = 2.35 (coefficient) × 2000 g (sample)</p>	<p>Put a sample to weigh on the balance.</p> <p>The weight of the sample is multiplied by the coefficient, and the result is displayed.</p>

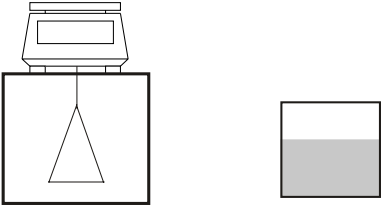

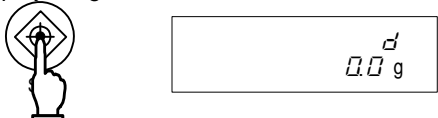
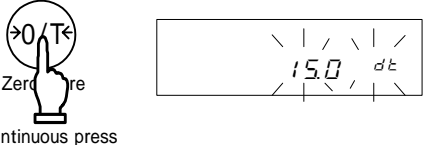
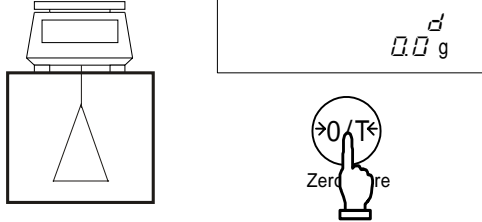
You can use the arrow keys to enter values instead of the **Zero/Tare** key or the **Function** key (Refer to page 22).

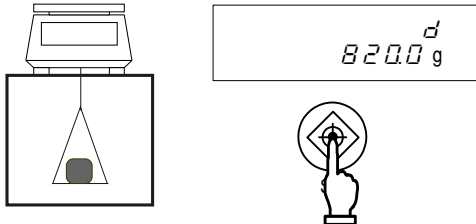
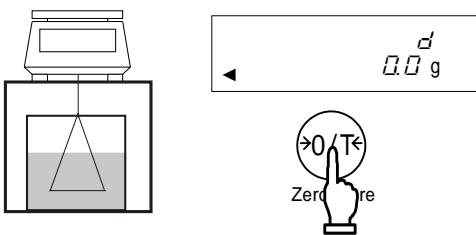
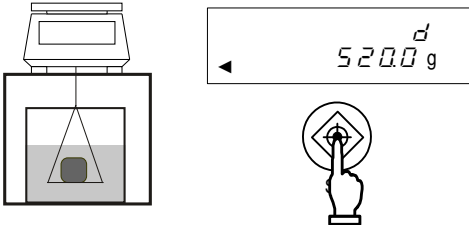
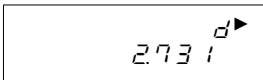
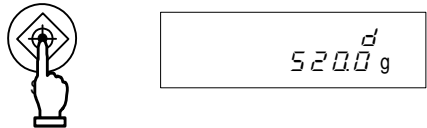
## 6.6 Gravimeter

When the balance is used as a gravimeter, specific gravity is measured by hanging a sample, thus requiring an optional hooking device. Any net or cage, hanging strings, or containers that may be required to hang the sample on the balance should be prepared by the user for the particular sample to be subjected to the measurement.

The measurement method complies with the requirements specified in JIS Z 8807, "Measuring Methods for Specific Gravity of Solid."

### 6.6.1 Measurement procedures for specific gravity

<p>1 Preparation for specific gravity measurements</p>  <p>Hanging measurement    Container with water (media)</p>	<p>Put the balance on a pedestal and attach the hanging hook. Hang a cage or the like that is used for holding the sample.</p> <p>As the sample weight is also measured in water, prepare a container with water whose size is enough to put the entire cage into the water.</p>
<p>2 Set the balance for gravimeter and select the media.</p> 	<p>Select the gravimeter function [ 1 SEt. 5 ] in Function 1.</p> <p>Select [ 0 : water ] or [ 1 : Other than water ] for the media [ 1.1 RE d. 0 ].</p>
<p>3 Display weight</p> 	<p>Press the <b>Set</b> key to show the weight.</p> <p>With gravimeter measurements, [d] will be indicated above the [g] symbol.</p>
<p>4. Measure the water temperature in the container and enter the value into the balance.</p>  <p>Continuous press</p>	<p>Press and hold the <b>Zero/Tare</b> key to enter the temperature setting mode and then enter the water temperature (for any media other than water, enter the specific gravity value of the media). (Refer to the descriptions on page 26 for the procedures to enter the value.)</p> <p>The value entered is held even after the power is turned off.</p>
<p>5 Reset the display to 0.</p> 	<p>Press the <b>Zero/Tare</b> key while only the cage is hung to reset the display to zero.</p>

<p>6. Measure the weight of the sample in the air.</p> 	<p>Measure the weight of the sample put in a cage or on a pan, in the air.</p> <p>After the weight display is stable, press the <b>Set</b> key to enter the weight.</p> <p>When the weight is entered, a [◀] is displayed in the bottom left of the display and a [ ] is shown in the top left part of the display.</p>
<p>7. Submerge only the hanging cage and set to zero.</p> 	<p>Before the measurement in water, submerge only the cage and set the indication to zero.</p> <p>This is for removing any residual error due to the hanging cage.</p> <p>* Pressing the <b>Print</b> key after entering the sample weight in the air will clear the weight value entered.</p>
<p>8. Measure the weight of the sample in water.</p> 	<p>Put the sample on the hanging cage and submerge the entire volume.</p> <p>Be careful that the cage does not touch the bottom of the container.</p> <p>After the weight display is stable, press the <b>Set</b> key to enter the weight.</p>
<p>9. The specific gravity value will be displayed.</p> 	<p>When the weight is entered, the resulting specific gravity value will be displayed.</p> <p>The [▶] symbol will be lit in the top right of the display while a specific gravity value is displayed.</p> <p>Specific gravity is has no unit.</p>
<p>10. Return to weight display by the <b>Set</b> key.</p> 	<p>When you press the <b>Set</b> key while a specific gravity value is displayed, the display returns to the weight display of the gravimeter function.</p>

### ☆ Hints ☆

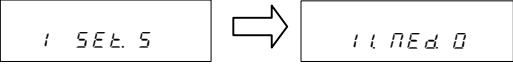
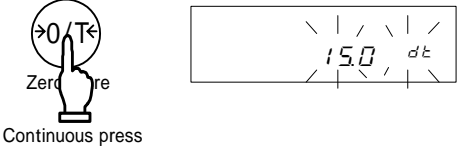
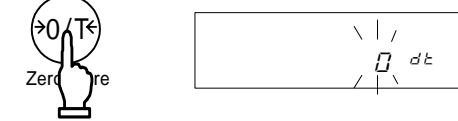
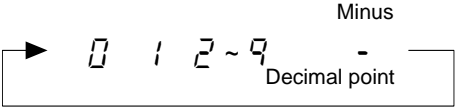

Measurements with smaller size samples may result in some variations in the measured specific gravity values. In particular, any sample whose specific gravity is larger tends to show more variation.

Use as large a sample as possible for obtaining more stable measurements.

## 6.6.2 Entering water temperature or the specific gravity of the media

This is for entering correction data for the liquid media to be used for the measurement. Enter water temperature if the media is water. Otherwise, enter the specific gravity value of the liquid. Use Function 1 for selecting the media as follows:

### Correction data input procedures

<p>1 Set the balance for gravimeter and select the media.</p> 	<p>Set Function 1 to select gravimeter and select either [ ] : water] or [ ] : Other than water] for the media [ ] : 11.17Ed. 0].</p>
<p>2 Call the mode to enter a correction value.</p> 	<p>Press and hold the <b>Zero/Tare</b> key to enter the temperature setting mode and then enter the water temperature (for any media other than water, enter the specific gravity value of the media). If the media is water, [d.t.] is displayed on the right side of the display. Otherwise, [d. ] will be indicated instead.</p>
<p>3 Press the <b>Zero/Tare</b> key.</p> 	<p>Press the <b>Zero/Tare</b> key. The [ ] indication blinks.</p>
<p>4. Enter the value.</p> 	<p>Use the <b>Zero/Tare</b> and <b>Function</b> keys for this entry. The <b>Zero/Tare</b> key is used for incrementing the figure and the <b>Function</b> key is used for moving the digit to the right. (This operation can also be made by using the arrow keys. Refer to page 22).</p>
<p>5 Save the value.</p> 	<p>Upon completing the entry, press the <b>Set</b> key to return to the weight display mode.</p>

### Caution

The acceptable range of numeric entry is specified as follows:

Media	Input Data	Range
Water	Water temperature	0.0 to 99.9°C
Any liquid other than water	Density	0.0001 to 9.9999 g/cm <sup>3</sup>

If any value outside the above specified range is entered, either the minimum or the maximum value within the range will be set.

The value set is held even after the power is turned off.

### 6.6.3 Gravimeter measurement data output

(1) Output before gravimeter measurement

With the gravimeter function, the operation is with the [5 I.O.C. 7] (output once when **Print** key is pressed after it is stabilized) mode, regardless of the setting by [5 I.O.C.].

(2) Output while specific gravity is displayed

[12 O.D.] is used for selecting the output contents.

The output method can also be selected by [13 R.O.].

(3) Output (print) format

When "specific gravity, weight, and actual water temperature or density of media" are all selected in [12 O.D.], .... If "specific gravity only" is selected, only the second line value will be output.

In addition, if any statistical calculation is made for the printer, the values in the second line will be printed being prefixed by a serial number.

Language selection between English and Japanese for printing is made by the GLP function [E 3 P.F.] of Function 1.

1) When water is selected

English

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15
1	D	E	N	S	I	T	Y		S	O	L	I	D		
2								x	x	.	x	x	x		
3	S	A	M	P	L	E		W	E	I	G	H	T		
4					x	x	x	x	x	.	x	x	x		g
5	T	E	M	P	E	R	A	T	U	R	E	N	O	W	
6					x	x	x	x	x	.	x				C
7															

Japanese (Katakana)

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15
1	コ	タ	イ	ヒ	シ	"	ユ	ウ							
2									x	x	.	x	x	x	
3	シ	"	ユ	ウ	リ	ヨ	ウ								
4					x	x	x	x	x	.	x	x	x		g
5	シ	"	ツ	ス	イ	オ	ン								
6					x	x	x	x	x	.	x				C
7															

2) When other than water is selected

English

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15
1	D	E	N	S	I	T	Y		S	O	L	I	D		
2								x	x	.	x	x	x		
3	S	A	M	P	L	E		W	E	I	G	H	T		
4					x	x	x	x	x	.	x	x	x		g
5	D	E	N	S	I	T	Y		M	E	D	.	L	I	Q
6					x	x	x	x	.	x	x	x			
7															

Japanese (Katakana)

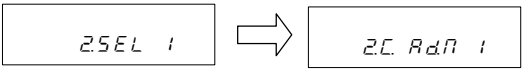

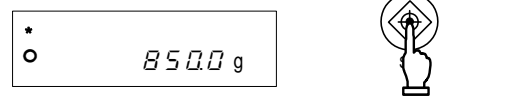
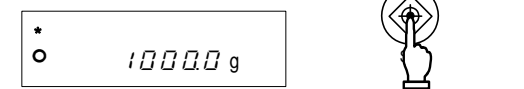
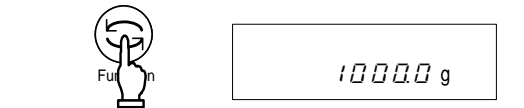
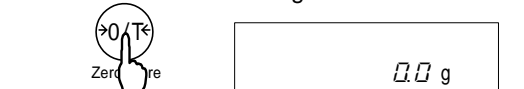
	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15
1	コ	タ	イ	ヒ	シ	"	ユ	ウ							
2									x	x	.	x	x	x	
3	シ	"	ユ	ウ	リ	ヨ	ウ								
4					x	x	x	x	x	.	x	x	x		g
5	ハ	"	イ	タ	イ	ミ	ツ	ト	"						
6					x	x	x	x	.	x	x	x	x		d s
7															

## 7 Addition Function

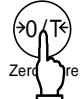
The addition function allows you to find the cumulative weight of samples put on the balance one after another. Except for the gravimeter function, all other functions, i.e., weighing machine, parts counting, percentage weighing, and unit converting functions support the addition function.

### 7.1 Cumulate Function

With the cumulate function, a cumulative weight can be obtained while samples are reloaded.

<p>1 Set Function 1.</p> 	<p>Set Function 1 from [2 SEL 1] to [2C Adn 1]. Then, press the <b>Set</b> key to return to the measurement mode.</p>
<p>2 Reset the display to 0.</p> 	<p>Press the <b>Zero/Tare</b> key for zero adjustment or for tare range.</p>
<p>3 Put samples on the balance, and obtain the weight.</p> 	<p>Put samples to weigh on the balance. When the asterisk (*) sign is displayed, you can put additional samples on. When [O] is displayed, press the <b>Set</b> key. The displayed value is added, and the cumulative value is displayed with the [Σ] sign for several seconds.</p>
<p>4 Put another sample on the balance.</p>	<p>Unload the balance, check that the display indicates 0, and put different samples on.</p>
<p>5 Read the weight.</p> 	<p>When [O] is displayed, press the <b>Set</b> key. After it is stabilized, the displayed value is added, and the cumulative weight is displayed with the [Σ] sign for several seconds.</p>
<p>6 Display the cumulative weight.</p> 	<p>Press the <b>Function</b> key twice. The sign [Σ] will be lit and the cumulative value will be shown. (The number of times the <b>Function</b> key is pressed depends on the type of the balance.)</p>
<p>7 Clear the cumulative weight.</p> 	<p>Press the <b>Zero/Tare</b> key with the cumulative weight displayed (step 4) to reset the display to 0.</p>

#### Obtain a cumulative weight without reloading the balance

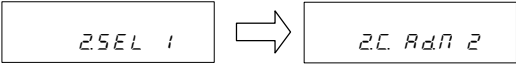
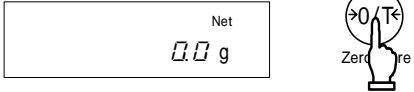
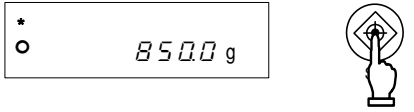
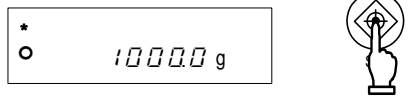
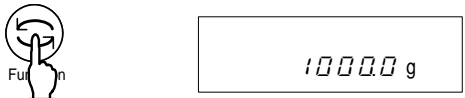
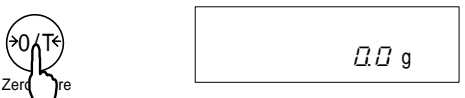
	<p>Press the <b>Zero/Tare</b> key in step 2 without unloading the balance. This will reset the display to 0. Put additional samples on, and press the <b>Set</b> key. Then, the weight is cumulated.</p>
---	--

#### Caution:

1. The addition function is available only when the display is once set to zero.
2. If [L - Err] is displayed when you press the **Set** key, it indicates that you put additional samples on twice, that you unloaded some samples, or that you pressed the key without adding samples.
3. You can add samples when an asterisk (\*) is displayed.
4. The [H L R] display with Function 1 can be used for turning ON/OFF the function to wait for stabilization upon additions.

## 7.2 Net Addition Function

The net addition function automatically performs tare range after an addition, so additions can be performed without having to replace samples.

<p>1 Set Function 1.</p> 	<p>Set Function 1 from [2.5EL 1] to [2.C. Adn 2]. Then, press the <b>Set</b> key to return to the measurement mode.</p>
<p>2 Reset the display to 0 by tare range.</p> 	<p>Press the <b>Zero/Tare</b> key for setting tare range.</p>
<p>3 Put samples on the balance, and obtain the weight.</p> 	<p>Put samples to weigh on the balance. When the asterisk (*) sign is displayed, you can put additional samples on. When [O] is displayed, press the <b>Set</b> key. The displayed value is added, and the cumulative value is displayed with the [Σ] sign for several seconds.</p>
<p>4 The balance automatically sets tare range.</p>	<p>The balance automatically sets the tare range.</p>
<p>5 Put additional samples on the balance to obtain the weight.</p> 	<p>When [O] is displayed after adding the sample, press the <b>Set</b> key. The displayed value is added, and the cumulative value is displayed with the [Σ] sign for several seconds.</p>
<p>6 Display the cumulative weight.</p> 	<p>Press the <b>Function</b> key twice. The sign [Σ] will be lit and the cumulative value will be shown. (The number of times the <b>Function</b> key is pressed depends on the type of the balance.)</p>
<p>7 Clear the cumulative weight.</p> 	<p>Press the <b>Zero/Tare</b> key with the cumulative weight displayed (step 6).</p> <p>* Any cumulative value will be cleared by a zero adjustment or tare range operation.</p>

### Caution

- When [E - Err] is displayed when you press the **Set** key, it indicates that you put additional samples on twice, that you unloaded some samples, or that you pressed the key without adding samples.
- You can add samples when an asterisk (\*) is displayed.

## 8 Limit Function

The limit function judges a weight based on limit values that you have stored on the balance.

The limit judgment result will be indicated by the position where the triangle [◀] is lit. When one or two points are set, the resulting display will be either HI (high), OK (proper), or LO (low). With a 3 or 4 point setting, 5 triangles [◀] are used to provide more detailed judgments.

### 8.1 Setting the Limit Function

Make the settings of the limit function in Function 1. While the limit function provides various setting items, carefully read Section 4.3: Additional Functions, on page 15, before actually performing the settings.

### 8.2 Judgment and Saving

A limit value can be judged in the following two ways. Select an appropriate one in Function 1.

- (1) Judge by absolute values..... Specify an upper weight limit or a lower weight limit directly.
- (2) Judge by deviation values..... Specify an upper limit or a lower limit with respect to a reference weight.

A limit value can be stored in the following two ways. They can be used in combination with either of the above judgment methods.

- (1) Put actual samples on the balance ..... Put an actual sample on the balance, and save the weight as a limit value.
- (2) Enter values ..... Enter a limit value using the keys.

- \* Once a limit value is entered, it is held after the balance is turned off.
- \* You can set separate limit values for the weighing machine, parts counting, percentage weighing, and unit converting functions.
- \* When you enter a limit value, it is treated differently between the two judgment methods. (Refer to Caution 9 on page 35.)

### 8.3 Display of Judgment Results

When one or two points are set, a [◀] is displayed next to one of HI, OK, and LO in the left of the display depending on the judgment result.

Judgment Result	When one point (lower limit) is set	When two points (lower and upper limits) are set
HI (high)	N/A	Weight > Upper limit
OK (proper)	Lower limit ≤ Weight	Lower limit ≤ Weight ≤ Upper limit
LO (low)	Lower limit > Weight	Lower limit > Weight

When 3 or 4 points are set, triangles [◀] are displayed in 4 or 5 levels according to the judgment result.

HI	◀	-----	Rank 5		
	◀	-----	Rank 4		
OK	◀	-----	Rank 3		
	◀	-----	Rank 2		
LO	◀	-----	Rank 1		


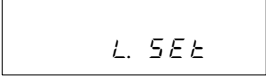
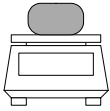

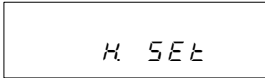


  

Judgment Result	3 or 4 points
Rank 5	Fourth point ≤ Weight
Rank 4	Third point ≤ Weight < Fourth point
Rank 3	Second point ≤ Weight < Third point
Rank 2	First point ≤ Weight < Second point
Rank 1	Weight < First point



## 8.4 Judge by Absolute Values

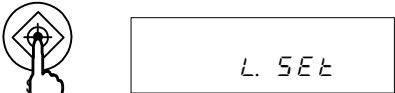

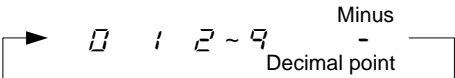


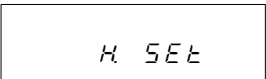

### 8.4.1 Set two limit values by putting actual samples on the balance

<p>1 Start the limit function.</p>  <p>Continuous press</p> 	<p>Press the <b>Set</b> key for a few seconds.</p> <p>When the display indicates [L. SET], release the key.</p> <p>The current lower limit is displayed next to [LO ◀].</p>
<p>2 Put the lower limit sample on the balance.</p> 	<p>Put the lower limit sample on the pan.</p>
<p>3 Save the lower limit.</p>  <p>Function</p>	<p>Press the <b>Function</b> key.</p> <p>When a lower limit is saved, the value is displayed for a while.</p> <p>* To set only one point, ignore the following steps.</p>
<p>4 Set an upper limit.</p> 	<p>The display then indicates [H. SET], and you can now set an upper limit.</p> <p>The current upper limit is displayed next to [HI ◀].</p>
<p>5 Put the upper limit sample on the balance.</p> 	<p>Put the upper limit sample on the pan.</p>
<p>6 Save the upper limit.</p>  <p>Function</p>	<p>Press the <b>Function</b> key.</p> <p>When an upper limit is saved, that value is displayed for a while. Then the balance goes back into measurement mode.</p>

\* For the case of three or four point setting, repeat the above steps 2 through 3. Limit values saved are displayed as [L 1 SET] - [L 3 SET] or [L 4 SET] rather than [L. SET] or [H. SET]. In addition, a triangle [◀] on the left side of the display will be lit at the same time for indicating the setting level.

HI	◀	-----	P L 4 SET	Fourth point
	◀	-----	P L 3 SET	Third point
OK				
	◀	-----	P L 2 SET	Second point
LO	◀	-----	P L 1 SET	First point

## 8.4.2 Set two limit values by entering values

<p>1 Start the limit function.</p>  <p>Continuous press</p>	<p>Press the <b>Set</b> key until the display indicates [L. SEt].</p> <p>The current lower limit is displayed next to [LO ◀].</p>
<p>2 Switch to the value input screen.</p>  <p>Zero/Tare</p>	<p>Press the <b>Zero/Tare</b> key.</p> <p>A (◻) blinks in the right end of the display.</p>
<p>3 Enter a number.</p>  <p>Minus - Decimal point</p>	<p>Press the <b>Zero/Tare</b> key to change the blinking number. Every time you press the key, the number changes.</p>
<p>4 Select digit places.</p>  <p>Function</p>	<p>Press the <b>Function</b> key, and the specified number shifts to the left, and you can enter the next digit.</p>
<p>5 Save the lower limit.</p> 	<p>Press the <b>Set</b> key.</p> <p>When a lower limit is saved, the value is displayed for a while.</p> <p>(To set only one point, ignore the following steps.)</p>
<p>6 Set an upper limit.</p> 	<p>The display then indicates [H. SEt], and you can now set an upper limit.</p> <p>The current upper limit is displayed next to [HI ◀].</p>
<p>7 Enter an upper limit.</p>	<p>Repeat the steps 2 to 5 to enter an upper limit.</p>
<p>8 Save the upper limit.</p> 	<p>Press the <b>Set</b> key. The upper limit is saved, and the balance goes back into measurement mode.</p>

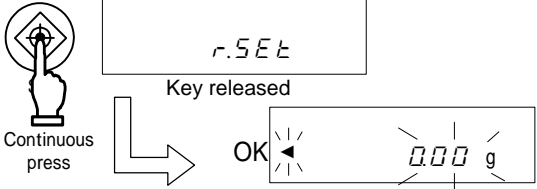
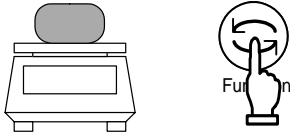
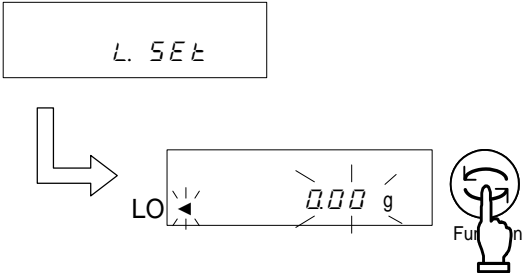
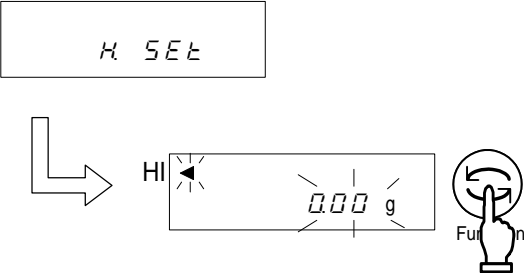
You can use the arrow keys to enter values instead of the **Zero/Tare** key or the **Function** key (Refer to page 23).

\* For the case of three or four point setting, repeat the above steps 2 through 5.

Limit values saved are displayed as [L 1 SEt] - [L 3 SEt] or [L 4 SEt] rather than [L. SEt] or [H. SEt]. In addition, a triangle [◀] on the left side of the display will be lit at the same time for indicating the setting level.

## 8.5 Judge by Deviation Values

### 8.5.1 Set two limit values by putting actual samples on the balance

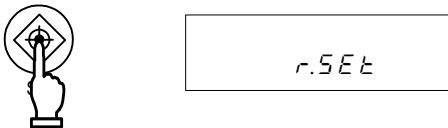
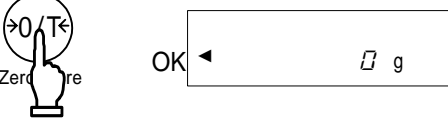
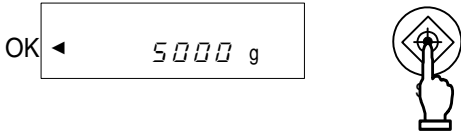
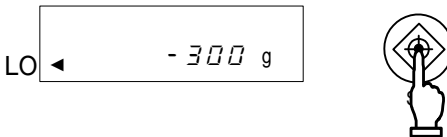
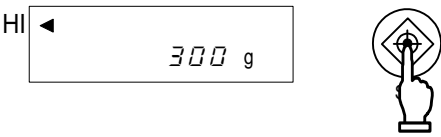
<p>1 Start the limit function.</p>  <p>Continuous press</p> <p>Key released</p> <p>OK</p>	<p>Press the <b>Set</b> key for a few seconds. When [r. SEt] is displayed, release the key.</p> <p>The current reference weight blinks in the display.</p>
<p>2 Save the reference weight.</p>  <p>Function</p>	<p>Put a sample that provides a reference weight on the balance, and press the <b>Function</b> key.</p> <p>When a reference weight is saved, that value is displayed for a while.</p>
<p>3 Save a lower limit.</p>  <p>L. SEt</p> <p>LO</p> <p>Function</p>	<p>The display first indicates [L. SEt.], and then the current lower limit blinks.</p> <p>Put a sample that provides a lower limit on the balance, and press the <b>Function</b> key. A lower limit is saved. Then go to the next step.</p>
<p>4 Save an upper limit.</p>  <p>H. SEt</p> <p>HI</p> <p>Function</p>	<p>The display first indicates [H. SEt], and then the current upper limit blinks.</p> <p>Put a sample that provides an upper limit on the balance, and press the <b>Function</b> key. The upper limit is saved, and the balance goes back into measurement mode.</p>

\* For the case of three or four point setting, repeat the above steps 3 and 4.

Limit values saved are displayed as [L 1 SEt] - [L 3 SEt] or [L 4 SEt] rather than [L. SEt] or [H. SEt]. In addition, a triangle [◀] on the left side of the display will be lit at the same time for indicating the setting level.

HI	◀	-----	PL 4 SEt	Fourth point
	◀	-----	PL 3 SEt	Third point
OK	◀	-----	r. SEt	Reference value
	◀	-----	PL 2 SEt	Second point
LO	◀	-----	PL 1 SEt	First point

## 8.5.2 Set two limit values by entering values

<p>1 Start the limit function.</p>  <p>Continuous press</p>	<p>Press the <b>Set</b> key for a few seconds. When [<i>r. S E t</i>] is displayed, release the key.</p> <p>The current reference weight blinks in the display.</p>
<p>2 Switch to the value input screen.</p> 	<p>Press the <b>Zero/Tare</b> key. A [<i>0</i>] blinks in the right end of the display.</p>
<p>3 Enter a reference weight.</p> 	<p>Perform the same procedures as given in steps 3 and 4 in Section 8.4.2, "Set two limit values by entering values" in page 32 to enter a reference weight.</p> <p>When you are entering a reference, the display indicates [<i>OK ◀</i>]. When you have entered a reference weight, press the <b>Set</b> key to save the value.</p>
<p>4 Save the lower limit.</p> 	<p>Repeat step 3 above to set a lower limit.</p> <div style="border: 1px solid black; padding: 5px;"> <p>To judge by deviation values, enter a value with a minus sign (-) to specify a difference from a reference weight (Refer to Caution 9 on page 35).</p> </div> <p>When you have entered the lower limit value, press the <b>Set</b> key to save the value. (To set only one point, ignore the following step and it returns to the measurement mode.)</p>
<p>5 Save the upper limit.</p> 	<p>Set an upper limit in the same way as in step 4 above.</p> <p>When you have entered an upper limit, press the <b>Set</b> key. The value is saved, and the balance goes back into the measurement mode.</p>

You can use the arrow keys to enter values instead of the **Zero/Tare** key or the **Function** key (Refer to page 22).

\* For the case of three or four point setting, repeat the above steps 2 through 5.  
Limit values saved are displayed as [*L 1 S E t*] - [*L 3 S E t*] or [*L 4 S E t*] rather than [*L S E t*] or [*H S E t*]. In addition, a triangle [*◀*] will be lit at the same time for indicating the setting level.

## Caution:

1. The initial limit values are all zero.
2. You can set separate limit values for the weighing machine, parts counting, percentage weighing, and unit converting functions. However, absolute values and deviation values are saved in the same storage area. So if you switch between absolute values and deviation values, limit values are deleted.
3. When the balance is not in measurement mode, you cannot set limit values by key operation. When a cumulative sum is being displayed, for example, you cannot set limit values.
4. Before you set limit values, adjust the zero-point or set a tare range as required (You can adjust the zero-point or set a tare range in sampling for parts counting or weight percent measuring).
5. You can check limit values by pressing the **Set** key.  
The lower limit is displayed after the display indicates [**L** **SET**], and the upper limit is displayed after [**H** **SET**].  
When you set three or four points, saved limit values are displayed as [**L 1 SET**] - [**L 3 SET**] or [**L 4 SET**] rather than [**L SET**] or [**H SET**].
6. When you made a mistake in an operation, press the **Print** key.  
As the operation will be cancelled, start the operation from the beginning.
7. If you press the **Function** key when a value is blinking in the display, the weight of the sample on the balance is set as an actual sample weight. If you press the **Zero/Tare** key here, you can now enter a limit value.
8. If the limit value entries are not lined up in the order of magnitude, five triangles [◀] will be lit regardless of the number of judgment points. Check the values and enter them again.
9. When you enter a limit value, it is treated differently between the two judgment methods.  
To judge by absolute values, enter limit values directly. On the other hand, to judge by deviation values, you have to specify a range with respect to a reference weight.

(Example)

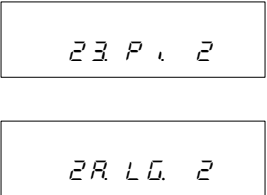
To set a lower limit of 970.0 g and an upper limit of 1050.0 g with respect to a reference weight of 1000.0 g, enter the limit values as shown below:

	Reference Weight	Lower Limit	Upper Limit
Absolute weight	1000.0 g	970.0 g	1050.0 g
Judgment by absolute values	1000.0 g	970.0 g	1050.0 g
Judgment by deviation values	1000.0 g	-30.0 g	50.0 g




## 8.6 Bar Graph for the 2-point Scale

You can set two points for the limit function, and display a bar graph to indicate the result in a range specified by the two points.

Except for the gravimeter function, all other functions, i.e., weighing machine, parts counting, percentage weighing, and unit converting functions can support this function.

<p>1 Set Function 1</p> 	<p>Set the 2-point scale [2] for Point Scale [2 3 P 1], and Bar graph for the 2-point scale [2] for How to Display Results [2 R LG] in Function 1.</p>
<p>2 Set upper and lower limits.</p>	<p>You can put actual samples on the balance or enter values to set upper and lower limits. You can also judge by absolute values or deviation values.</p>

A bar graph is displayed in one of the following ways:


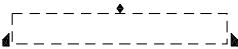
Bar graph	Weight range	Description
	$\text{Lower limit} > \text{Weight}$	A bar graph is not displayed.
	$\text{Lower limit} \leq \text{Weight} \leq \text{Upper limit}$	A bar graph to indicate the judgment result is displayed.
	$\text{Weight} > \text{Upper limit}$	The whole bar graph is displayed.

### Caution:

- When the lower limit is the same as the upper limit, no bar graph will be displayed.
- When a 2-point bar graph is being displayed, a normal bar graph to indicate the rate of a weight to the weighing capacity is not displayed.
- When a 2-point bar graph is being displayed, the limit function does not work.

### ☆ Hints ☆

Different bar graph frames are used depending on whether the balance is in the measurement mode or in the 2-point bar graph display mode, as shown below:

Bar graph frame	Balance condition
	Measurement mode
	2-point bar graph

## 9 Calibration and Span Test for the Balance

### 9.1 Span Adjustment

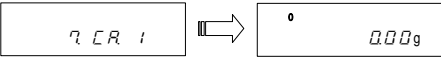
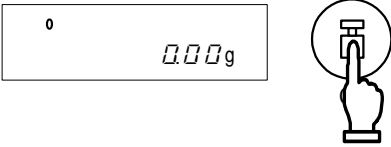
An electronic balance, which is influenced by the acceleration of gravity, indicates different values depending on the location it is used. For this reason, you should calibrate your balance every time you relocate it. You should also calibrate it after a long time of no use or when it does not indicate correct values.

To calibrate a balance is called span adjustment, which is required for highly accurate measurements.

- \* You must wait for at least 30 minutes after turning the power on, before performing span adjustment.

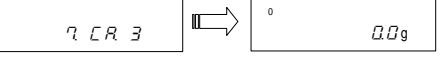
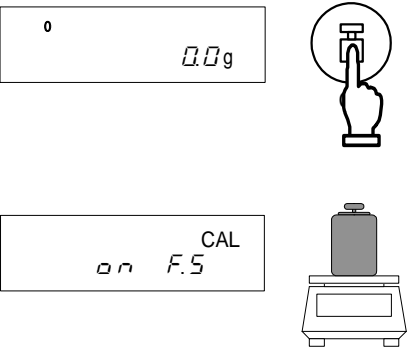
#### 9.1.1 Span Adjustment with Built-in Weight (Auto-span Adjustment)

\* Only supported by the HJR series units

<p>1 Set Function 1.</p> 	<p>Select span adjustment with built-in weight (auto-span adjustment) in Function 1 [FUNCTION 1]. Then, press the <b>Set</b> key to return to the weight display.</p>
<p>2 Start span adjustment.</p> 	<p>Press the <b>Cal</b> key.</p> <p>The built-in weight is enabled, and automatically starts span adjustment. The display indicates [Auto CAL], [WRT]*, [CHK 0], [CHK F.S.], [BUSY], and [END] in sequence. When span adjustment is completed, the balance goes back to the normal weight display.</p>

- \* While [WRT] is displayed, the balance is preparing to read the zero point. The remaining time is indicated by the bar graph display. (Displayed with HJR-2200E, HJR-3200E, HJR-4200E)


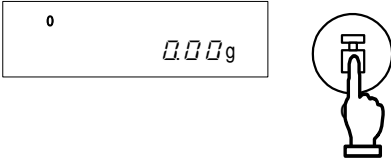
#### 9.1.2 Span Adjustment with External Weight

<p>1 Set Function 1.</p> 	<p>Select span adjustment with built-in weight (auto-span adjustment) in Function 1 [FUNCTION 3]. Then, press the <b>Set</b> key to return to the weight display.</p>
<p>2 Start span adjustment.</p> 	<p>Press the <b>Cal</b> key.</p> <p>The display first indicates [CAL ERR], and then [0.0 0], and the balance starts adjusting the zero-point. (If weight error is used, the corresponding weight error entered will be indicated for about 2 seconds.)</p> <p>When the display indicates [0.0 F.S.], put a weight on the pan. The balance starts adjusting the weighing capacity.</p> <p>The display indicates [BUSY] and then [END] in sequence. When span adjustment is completed, the balance goes back to the normal display.</p>

## 9.2 Span Test


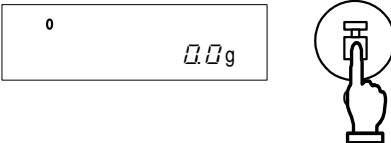
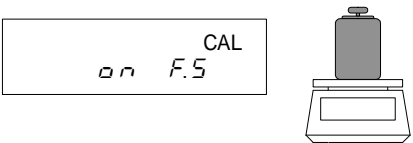
The purpose of the span test is to check for any offset in the span of the balance with respect to a reference weight. Performing this function will not calibrate the balance.

### 9.2.1 Span Test with Built-in Weight \* Only supported by the HJR series

<p>1 Set Function 1.</p> 	<p>Set Function 1 with nothing put on the pan. Select span test with built-in weight in Function 1 [FUNCTION 2]. Then, press the <b>Set</b> key to return to the weight display.</p>
<p>2 Start the span test.</p> 	<p>Press the <b>Cal</b> key.</p> <p>The built-in weight is enabled, and automatically starts a span test. The display indicates [E. INT], [UR 1E]*, [E. 0], [E. F.S], [d IFF], and an error in sequence. When you press a key, the balance goes back into the measurement mode.</p>

\* While [UR 1E] is displayed, the balance is preparing to read the zero point. The remaining time is indicated by the bar graph display. (Displayed with HJR-2200E, HJR-3200E, HJR-4200E)

### 9.2.2 Span Test with External Weight

<p>1 Set Function 1.</p> 	<p>Set Function 1 with nothing put on the pan. Select span test with external weight in Function 1 [FUNCTION 4]. Then, press the <b>Set</b> key to return to the weight display.</p>
<p>2 Start a span test.</p>  	<p>Press the <b>Cal</b> key.</p> <p>The display first indicates [E. EHE], and then [0.0 0], and the balance starts testing the zero-point. (If weight error is used, the corresponding weight error entered will be indicated for about 2 seconds.)</p> <p>When the display indicates [0.0 F.S], put a weight on the pan. The balance starts a weighing capacity point test.</p> <p>The display first indicates [d IFF] and then an error. When you press a key, the balance goes back into the measurement mode.</p>

\* An error displayed after [d IFF] is calculated as shown below:


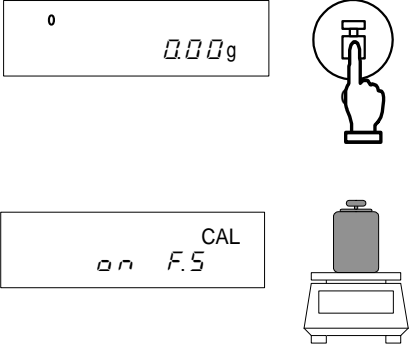
$$\text{Error} = \text{True value} - \text{Current weight}$$

A positive value displayed as an error means that the balance indicates a weight lighter than the actual weight.



### 9.3 Calibration of Built-in Weight \* Only for the HJR series

This is a function to calibrate the built-in weight with an external weight.

<p>1 Set Function 2.</p> 	<p>Set Function 1 with nothing put on the pan.          Enable the calibration of built-in weight function in Function 2 [3 r.ER 1].          Then, press the Set key to return to the weight display.</p>
<p>2 Start span adjustment.</p> 	<p>Press the <b>Cal</b> key. The display is changed [r EF. CAL].          Press the <b>Function</b> key while holding down the <b>Zero/Tare</b> key, and release the keys at the same time. The display changes from [UR It] to blinking [0 n 0], and then indicates [0 n F.5].            When the display indicates [0 n F.5], put a weight on the pan carefully. The display blinks [0 n F.5], and then indicates [0 n 0].            When the display is changed from [UR It] to [0 n 0], unload the weight. When calibration is completed, the display indicates [Er d], and the balance goes back to measurement mode.</p>

\* While [UR It] is displayed, the balance is preparing to read the zero point. The remaining time is indicated by the bar graph display. (Displayed with HJR-2200E, HJR-3200E, HJR-4200E)

#### ☆ Hints ☆

1. When you press a **key other than the Function key** during span adjustment or a span test, the display indicates [St OP], and the span adjustment or test is cancelled. Then the balance goes back to the normal display.
2. To perform span adjustment or a span test with an external weight, use a weight for calibration that weighs 50% of the weighing capacity or heavier. To calibrate more accurately, use a weight that is close to the weighing capacity.
3. When an error is found in span adjustment or in a span test, the following messages are displayed:

- (1) [1 - Err]: You used a weight weighing less than 50% of the weighing capacity for span adjustment.
- (2) [2 - Err]: An error over 1.0% was detected in span adjustment with an external weight, or the balance failed.
- (3) [3 - Err]: Automatic span adjustment was performed with an object put on the pan.
- (4) [4 - Err]: An error over 1.0% was detected in automatic span adjustment, or the balance failed.
- (5) [R - Err]: Abnormal operation occurred in the internal driving unit during automatic span adjustment.

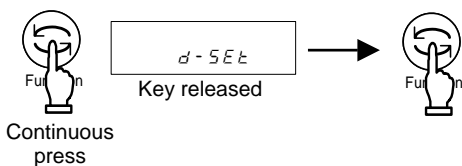
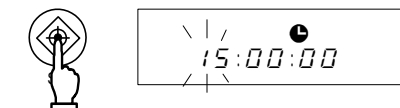

(Pressing a key when an error message is being displayed moves you back to measurement mode).

\* When one of these messages is displayed, calibration will not be performed. Check that you are using a proper weight, and start span adjustment or a span test again. If span adjustment with a proper weight results in the same error message, contact our Marketing Division or Technical Service Division.

# 10 Date and Time Setup

## 10.1 Time Setup

The time is displayed with a clock [🕒] symbol. Set the time in Hour-Minute-Second format on a 24-hour basis.

<p>1 Display the time.</p>  <p>Continuous press</p>	<p>Press the <b>Function</b> key for a few seconds. When the display is changed from [F u n c] to [d - 5 E t], release the key.</p> <p>Press the <b>Function</b> key one time. The display automatically indicates the time after indicating [t i m e].</p>
<p>2 Set the time.</p> 	<p>Press the <b>Set</b> key. You can change a blinking digit.</p> <p>Use the <b>Zero/Tare</b> key to change the blinking digit. Use the <b>Function</b> key to blink the next digit to the right. If you press the <b>Function</b> key when the rightmost digit is blinking, the leftmost digit blinks.</p>
<p>3 Save the time.</p> 	<p>Press the <b>Set</b> key to save the settings. The changes are saved, and a date is displayed.</p>

When you press the **Print** key before saving a time, the original time is displayed.

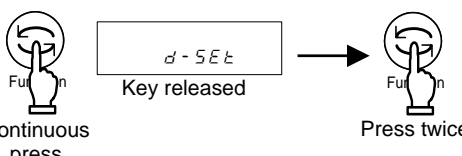
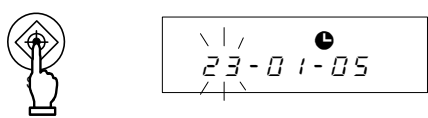

### 30-sec adjustment function

Press the **Zero/Tare** key when the time is being displayed.

29 or fewer seconds are rounded down to the nearest minute, and 30 or more seconds are rounded up to the nearest minute.

## 10.2 Date Setup

A date is displayed with a clock [🕒] symbol. The last two digits of the year are displayed. You can change the date format in Function 1.

<p>1 Display the date.</p>  <p>Continuous press</p>	<p>Follow step 1 in 10.1 to change the display to [d - 5 E t]. (Go straight to step 2 if you set a date after completing steps in 10.1.)</p> <p>Press the <b>Function</b> key twice. The display automatically indicates a date after indicating [d R t E].</p>
<p>2 Set a date.</p> 	<p>Press the <b>Set</b> key. You can change a blinking digit.</p> <p>Use the <b>Zero/Tare</b> key to change the blinking digit. Use the <b>Function</b> key to blink the next digit to the right. If you press the <b>Function</b> key when the rightmost digit is blinking, the leftmost digit blinks.</p>
<p>3 Save the date.</p> 	<p>Press the <b>Set</b> key to save the settings. The changes are saved, and the balance goes back into measurement mode.</p>

When you press the **Print** key before saving a date, the original date is displayed.

# 11 Miscellaneous Functions

---

## 11.1 Auto Sleep Function

This is a function to turn off the display when the balance is left to stand in measurement mode for about 3 minutes. This function can only be used when the balance is operated on the AC adapter. To use the auto sleep function, enable it in Function 1 [*R* *R.S.* *I*].

When the auto sleep function is enabled, the Sleep lamp (LED) lights up.

To exit the sleep mode, touch the pan or press a key.

### **Caution:**

The auto sleep function does not work under the following conditions:

1. You are setting Function 1, Function 2, time or date, or the interval function.
2. An object is put on the pan, and the display is not stable.
3. The balance is operated on batteries.

## 11.2 Auto Power Off Function

This is a function to turn off the balance when the balance is left to stand in measurement mode for about 3 minutes. This function can only be used when the balance is operated on batteries. This function saves battery capacity.

To use this function, enable it in Function 1 [*Q* *R.P.* *I*].

### **Caution:**

The auto power off function does not work under the following conditions:

1. You are setting Function 1, Function 2, time or date, or the interval function.
2. An object is put on the pan, and the display is not stable.
3. The balance is operated on the AC adapter.

## 11.3 Set Unit Function

This is a function to allow you to select two units (unit A and unit B) and switch between the units.

Which unit to display can be set in Function 1. You can also set the same unit for unit A and unit B.

Pressing the Function key switches the unit between unit A and unit B.

### **Caution:**

To use unit B, select the weighing machine function in Function 1 [*I* *SEt* *I*]. Unit A can be used in all the measurement types of the balance.

## 11.4 Minimum Readability Setting Function

This is a function to change the minimum readability. Set a minimum readability in Function 1. When you set a rougher or larger minimum readability, weighing is performed faster. The minimum readability depends on the unit.

You can also set the same unit for unit A and unit B in 11.3, and switch the minimum readability between unit A and unit B.

\* The least significant digit cannot be set higher than the 10's digit.

## 11.5 Advice CAL

Advice CAL blinks the CAL sign in the display to inform you when span adjustment is to be made. The timing to perform span adjustment depends on changes in ambient temperature and used hours. To use this function, set Function 1 to [*d. R d. L. 1*]. When the CAL sign starts to blink, perform span adjustment as soon as possible.

## 11.6 Date Display

This is a function to set the date format to be displayed in the balance or output to a printer. Set date format in Function 1 appropriately as shown below:

Date Display	<i>F. d R t E</i>	<i>1</i>	Output in Year-Month-Day format.
		<i>2</i>	Output in Month-Day-Year format.
		<i>☆ 3</i>	Output in Day-Month-Year format.

## 11.7 Time Stamp Output

This is a function to output measurement data with the current time. You should set the time before this function is used.

To use this function, enable it in Function 1 [*L. t. d. 1*].

## 11.8 Save Tare Weight Function

This is a function to use the latest tare weight when the **On/Off** key is pressed to turn on the balance. Use this function when you turn the balance on or off with a sample and tare put on the pan.

To use this function, enable it in Function 1 [*L. t R r E 1*]. The saved tare weight is updated every time a tare range is set.

### **Caution:**

Leaving the balance loaded with a substance and tare for a long period may result in a larger error in weighing. You should set the tare range regularly.

## 11.9 Direct Start Function




This is a function to turn on the balance automatically when it is connected to the AC power. You can use this function when the balance is used in conjunction with other devices. To use this function, enable it in Function 1 [*L. d S t. 1*].

## 11.10 Interval Output Function

This function outputs data at regular intervals. You can also output data time stamped. Set an interval in Hour-Minute-Second format. To use the interval output function, enable it in Function 1 to [S t a.c. R] or [S t a.c. b].

\* If it is desired to output time stamp data at the same time, refer to Section 11.7.

### 11.10.1 Set interval output

<p>1 Call the interval function.</p>  <p>Continuous press</p>	<p>Press the <b>Set</b> key until the display indicates [Int. URL]. The display indicates an interval with the leftmost digit blinking.</p>
<p>2 Set an interval.</p> 	<p>Press the <b>Zero/Tare</b> key to change the figure at the blinking digit. Use the <b>Function</b> key to move the blinking position to the right. If you press the <b>Function</b> key while the rightmost digit is blinking, the leftmost digit blinks. (You can use the arrow keys instead. Refer to page 22.)</p>
<p>3 Save the settings, and go back to measurement mode.</p> 	<p>Press the <b>Set</b> key to return to measurement mode.</p> <p>When you press the <b>Print</b> key before saving an interval, the original interval is displayed. If you press any key other than the <b>Set</b> key here, the balance goes back to the measurement mode, but the settings are not saved.</p>

### 11.10.2 Start interval output

Press the **Print** key. The display indicates [S t R r t], and starts interval output.

A [●] blinks during interval output. A [☞] lights up when data is output.

### 11.10.3 Stop interval output

Press the **Print** key. The display indicates [E n d], and the balance goes back to measurement mode.

#### Caution:

1. Note that data may not be output at exact intervals because output operation is interrupted while Function 1 or Function 2 is displayed.
2. If [E - E r r ] is displayed, the interval is set to zero. In this case, you should set an interval again.


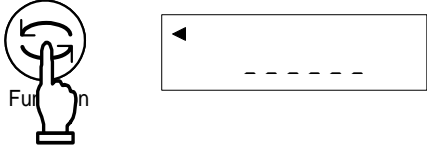
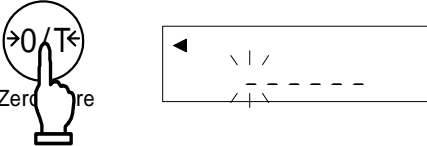
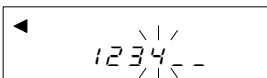

## 11.11 ID Number Entry

An ID number is used when data is printed in compliance with ISO/GLP/GMP. Set an ID number when you print data.

When an ID number is displayed, two triangles, [◀] and [ ] are displayed in the upper left part of the display.

You can use up to 6 digits in an ID number. You can use [0-9], [A-F], and [ - ], which are displayed in this order. An underscore [ \_ ] indicates a blank space.

### ID Number setting procedures


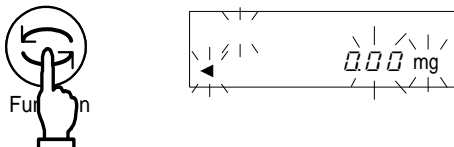
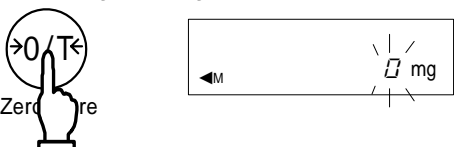
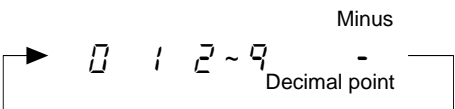

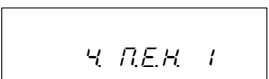
<p>1 Set Function 2.</p> 	<p>Enable Setup of ID No. in Function 2 [ i id i ].</p>
<p>2 Display an ID number.</p> 	<p>Press the <b>Function</b> key. An ID number is displayed.</p> <p>If you want to check the ID number, press the <b>Function</b> key again. The display is changed to [2 a.N.P. 0].</p>
<p>3 Start to enter an ID number.</p> 	<p>Press the <b>Zero/Tare</b> key. You can now enter an ID number.</p> <p>The leftmost digit blinks first.</p>
<p>4 Enter an ID number.</p> 	<p>Press the <b>Zero/Tare</b> key to change the blinking value, and the <b>Function</b> key to blink the next digit to the right. If you press the <b>Function</b> key when the rightmost digit is blinking, the leftmost digit blinks. (You can use the arrow keys instead. Refer to page 22.)</p>
<p>5 Save the ID number.</p> 	<p>Lastly, press the <b>Set</b> key. The setting is saved, and the display is changed to [2 a.N.P. 0].</p> <p>When you press the <b>Print</b> key before saving an ID number, the original ID number is displayed.</p>

## 11.12 Entry of Weight Error

By entering the weight error of an external reference weight to be used in span adjustment or span test, a more accurate calibration can be performed. Enter the weight error obtained from the following equation: (Weight errors are entered in the unit of mg.)

$$\text{Weight Error} = \text{Actual Weight} - \text{Nominal Value}$$

Example: For a 1000 g weight, a weight error may be  $= 1000.00012 - 1000$   
 $= 0.00012 = +0.12 \text{ mg}$

<p>1 Set Function 2.</p> 	<p>Set Function 2 to [2 a.N.P. 1].</p>
<p>2 Display the weight error.</p> 	<p>Press the <b>Function</b> key. The display shows a weight error.</p> <p>If you are only to check the weight error value, press the <b>Function</b> key again. The display is changed to [3 r.C.R. 0] or [4 N.E.H. 0] for the next item.</p>
<p>3 Start entering the weight error.</p> 	<p>Press the <b>Zero/Tare</b> key. You can now enter the weight error.</p> <p>The leftmost digit blinks first.</p>
<p>4 Entry of weight error.</p> 	<p>Press the <b>Zero/Tare</b> key to change the blinking value, and the <b>Function</b> key to blink the next digit to the right. If you press the <b>Function</b> key when the rightmost digit is blinking, the leftmost digit blinks. (You can use the arrow keys instead. Refer to page 22.)</p>
<p>5 Save the weight error data.</p> 	<p>Lastly, press the <b>Set</b> key. The setting is saved, and the display is changed to [3 r.C.R. 0] or [4 N.E.H. 0] for the next item.</p> <p>If you press the <b>Print</b> key during the entry, the display returns to the previous weight error value.</p>
<p>6 Use the weight error.</p> 	<p>For using the weight error entered, set to [4 N.E.H. 1].</p>

### ☆ Hints ☆

If more than one weight is used in combination, enter the total error of the weights to be used.

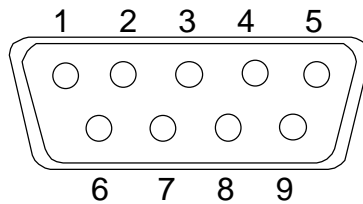
The range of the weight error should be within  $\pm 100.00 \text{ mg}$ . If any value outside the range is entered, it will result in the display of [r - E r r].

## 12 Input/Output Functions

### 12.1 RS232C Output

#### 12.1.1 Connector pin numbers and functions

Pin No.	Signal Name	Input/Output	Function & Remarks
1	-	-	-
2	RXD	Input	Receiving data
3	TXD	Output	Transmitting data
4	DTR	Output	HIGH (When the balance is powered ON)
5	GND	-	Signal ground
6	-	-	-
7	-	-	-
8	-	-	-
9	-	-	-



D-SUB9P Male Connector: Rear Panel

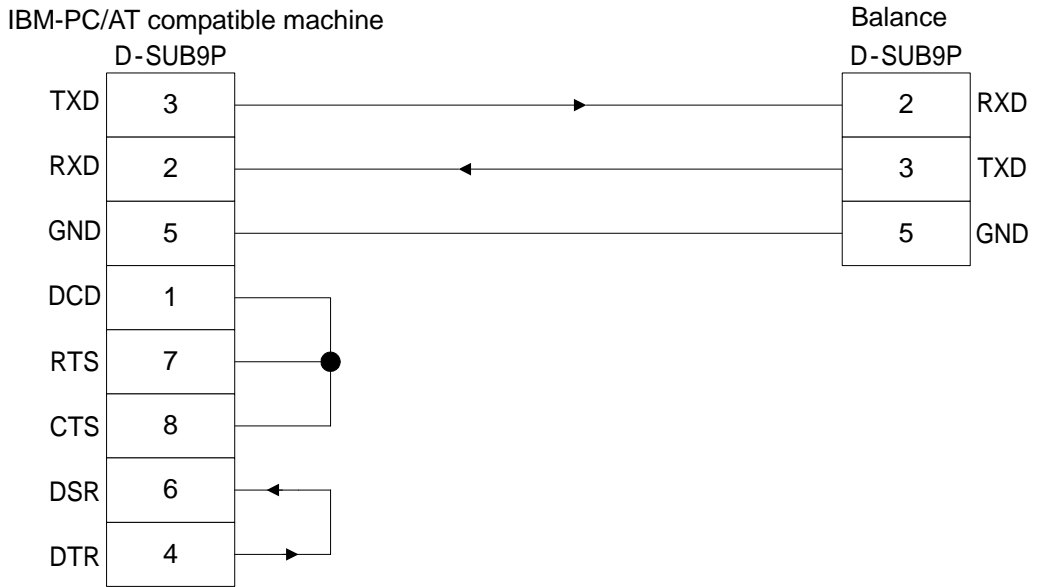
#### **Caution:**

Be sure to disconnect the AC adapter from the power outlet before plugging or unplugging the connector.

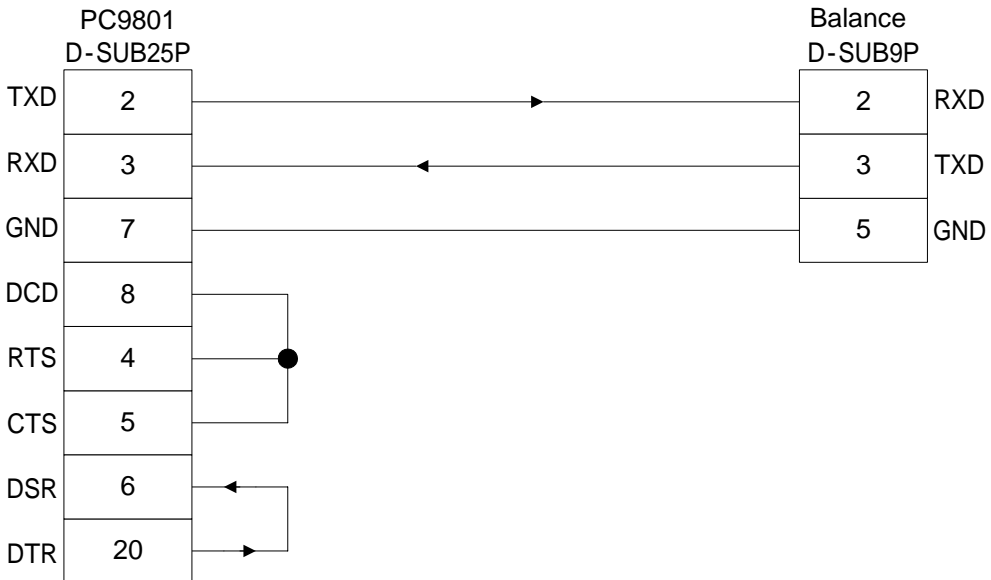


## 12.1.2 Connecting the balance to a PC

### ■■■ Sample connection with an IBM-PC/AT compatible machine ■■■

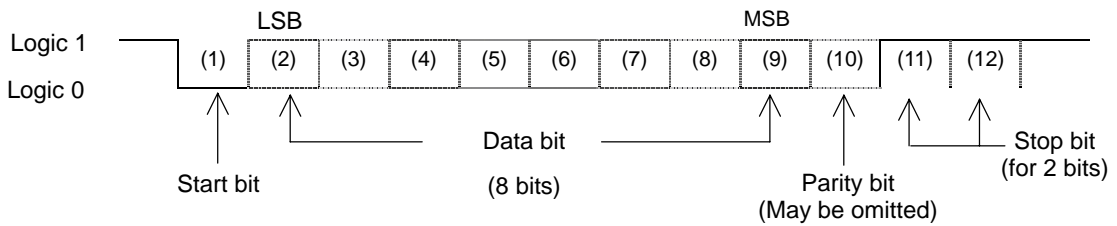


### ■■■ Sample wire connection with PC9801 ■■■



### 12.1.3 Interface specifications

- (1) Transmission system    Serial transmission, Start-stop synchronization
- (2) Transmission rate      1200/2400/4800/9600/19200 bps
- (3) Transmission codes    ASCII codes (8/7 bits)
- (4) Signal level            Compliant with EIA RS-232C.  
HIGH level (data logic 0) +5 to +15 V  
Low level (data logic 1) -5 to -15 V
- (5) Bit configuration      Start bit            1 bit  
Data bits            8/7 bits  
("7 bit" can be specified only for the extended 7-digit numeric format.)  
Parity bit           0/1 bit  
Stop bits            2/1 bit  
("1 bit" can be specified only for the extended 7-digit numeric format.)
- (6) Parity bit            None/Odd/Even



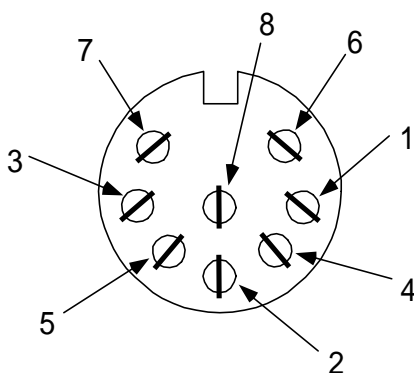
## 12.2 Output to Peripherals

Our standard peripheral units can be connected to the balance.\*<sup>1</sup>

These peripheral units include: CSP-160, CSP-240

### 12.2.1 Connector pin numbers and functions

Pin No.	Signal Name	Input/Output	Function & Remarks
1	EXT.TARE	Input	Tare setting from an external device* <sup>2</sup>
2	-	-	-
3	-	-	-
4	TXD	Output	Transmitting data
5	GND	-	Signal ground
6	-	-	-
7	-	-	-
8	-	-	-



DIN8P connector for output to peripheral devices

- \*1 To connect a peripheral unit to the balance, use the connection cable supplied with the unit.
- \*2 You can set a tare range or adjust the zero-point. from an external device by connecting a contact or a transistor switch between pin 1 and pin 5. In this case, allow at least 400 ms for connection (ON) time (Maximum voltage when the balance is turned OFF: 15 V, sink current when it is turned ON: 20 mA).

### **Caution:**

Be sure to disconnect the AC adapter from the power outlet before plugging or unplugging the connector.

## 12.3 Type of Communication Texts

This interface function uses the following three types of communication texts:

- (1) Output data            Data, such as weight values, that is output from the balance to an external unit
- (2) Input commands        Commands to control the balance from an external unit
- (3) Response                Response that is output from the balance to an input command

### Caution:

1. All of these communication texts can be used for RS232C output.
2. Only output data (1) can be used for a peripheral output unit.

## 12.4 Output Data

You can select one from the following three formats in Function 1:

### 12.4.1 Data format

#### (1) 6-digit numeric format

Consists of 14 characters including terminators (CR=0DH, LF=0AH).

1	2	3	4	5	6	7	8	9	10	11	12	13	14
P1	D1	D2	D3	D4	D5	D6	D7	U1	U2	S1	S2	CR	LF

#### (2) 7-digit numeric format

Consists of 15 characters including terminators (CR=0DH, LF=0AH). A parity bit can be appended.

1	2	3	4	5	6	7	8	9	10	11	12	13	14	15
P1	D1	D2	D3	D4	D5	D6	D7	D8	U1	U2	S1	S2	CR	LF

#### (3) Extended 7-digit numeric format

Consists of 15 characters including terminators (CR, LF). A parity bit can be appended.

This is an extended version of the 7-digit numeric format. The extended version is different from the 7-digit numeric format in that:

- the data length can be 7 bits rather than 8 bits, and
- the stop bit may be set to 1 bit rather than 2 bits.

The extended version is the same as the normal 7-digit format in other bits. However, if you select Japanese (Katakana) for Printed Language [*E 3 P.F.*] in Function 1, the data length is forced into 8 bits.

### 12.4.2 Polarity (P1: 1 character)

P1	Code	Description
+	2BH	When data is 0 or positive.
-	2DH	When data is negative.

### 12.4.3 Numeric data

**6-digit numeric format: (D1-D7: 7 characters)**

**7-digit numeric format: (D1-D8: 8 characters)**

D1 to D7(D8)	Code	Description
0 – 9	30H to 39H	Digits 0 – 9
•	2EH	* Decimal point (floating) * Omitted when numeric data does not contain decimal places. In this case, space is output to the least significant digit.
SP (space)	20H	* A space heading a numeric value * When numeric data does not contain decimal places, space rather than a decimal point is output to the least significant digit <sup>*1</sup> .

\*1 A numeric value is headed by 0 (30H) by factory default setting. You can change the function setting so that it is headed by space (20H) rather than 0.

### 12.4.4 Unit (U1, U2: 2 characters)

\* Codes are all in ASCII codes.

U1	U2	Code		Meaning	Balance indicators
(SP)	G	20H	47H	gram	<b>g</b>
K	G	4BH	47H	Kilogram	<b>kg</b>
C	T	43H	54H	carat	<b>ct</b>
O	Z	4FH	5AH	ounce	<b>oz</b>
L	B	4CH	42H	pound	<b>lb</b>
O	T	4FH	54H	troy ounce	<b>oz t</b>
D	W	44H	57H	pennyweight	<b>dwt</b>
G	R	4BH	52H	Grain	<b>GN</b>
T	L	54H	4CH	tael (Hong Kong)	<b>t</b>
T	L	54H	4CH	tael (Singapore, Malaysia)	<b>t</b> Top right ▶
T	L	54H	4CH	tael (Taiwan)	<b>t</b> Middle right ▶
M	O	4DH	4FH	momme	<b>mom</b>
t	o	74H	6FH	tola	<b>to</b>
P	C	50H	43H	Pieces	<b>Pcs</b>
(SP)	%	20H	25H	Percentage	<b>%</b>
(SP)	#	20H	23H	Computation results, numbers, etc.	<b>#</b>

#### 12.4.5 Judgment result when the limit function is enabled (S1: 1 character)

S1	Code	Description	
L	4CH	Too little (LO)	1- or 2-point scale
G	47H	Proper (OK)	
H	48H	Too much (HI)	
1	31H	Rank 1	3- or 4-point scale
2	32H	Rank 2	
3	33H	Rank 3	
4	34H	Rank 4	
5	35H	Rank 5	
T	54H	Cumulative value	Data type
U	55H	Unit weight	
(SP)	20H	No judgment result or no data type specified	
d	64H	Gross	

#### 12.4.6 Status (S2:1 character)

S2	Code	Description
S	53H	Data stable *
U	55H	Data unstable *
E	45H	Data error (Indicates that data other than S2 is invalid.) ( [a - E r r], [u - E r r] )
(SP)	20H	No status specified

\* Data that is independent of whether the weighing condition is stable or not, such as cumulative values and unit weights, is independent of whether S2 is S or U when it is output.

#### 12.4.7 Output data other than measurement data

The following data is accompanied by our proprietary printer control commands prior to and after the data when output.

##### (1) Interval output

When each interval output is started or stopped, a header and a footer are output.

##### Header

"-" is output for 15 characters.

1	2	3	4	5	6	7	8	9	10	11	12	13	14	15
-	-	-	-	-	-	-	-	-	-	-	-	-	-	-

##### Footer

Two linefeeds are inserted.

##### (2) Time output

When the time stamp function is enabled, the time is output one line above the output data.

1	2	3	4	5	6	7	8
h	h	:	m	m	:	s	s

\* hh: hour (00-23), mm: minute (00-59), ss: second (00-59)

## 12.5 Input Commands

The following 7 input commands are supported:

- |                                  |                                 |
|----------------------------------|---------------------------------|
| (1) Tare range command           | (2) Set output control command  |
| (3) Set measurement mode command | (4) Request date output command |
| (5) Request time output command  | (6) Set interval command        |
| (7) Span adjustment/test command |                                 |

### 12.5.1 Procedure for transmission

- (1) An input command is sent from an external device.

The full-duplex transmission system allows you to send an input command at any time independently of the data transmit timing of the balance.

- (2) Upon successful completion of an input command, the balance will send out either a normal completion response or the data requested by the command. If the operation has not resulted in normal completion, or if the received input command is invalid (in error), the balance will transmit an error response.

When the balance is in normal display mode, it usually sends a response to an input command within one second after it is received.

However, if the balance receives a tare range command when the Set Tare function is set so that the display is reset to 0 when the balance is stabilized after the Zero/Tare key is pressed, or if it receives an input command that takes a long time, it sends a response after the command is completely processed.

If the balance receives an input command when you are setting a function, when it is under span adjustment, or it is busy for other reasons, it sends a response after that operation is completed.

- 1) Commands to which a response is sent back after a specified operation is executed:
  - Tare range command, span adjustment/test command
- 2) Commands to which a response is sent back immediately or in 1 second after they are received:
  - Commands other than those specified in 1) above
- 3) Once you have sent an input command from an external device, do not send another input command to the balance until the external device receives a response to the command from the balance.

## 12.5.2 Response

You can select the response format of either the A00/Exx format or the ACK/NAK format in Function 1.

### (1) A00/Exx format

Consists of 5 characters including terminators (CR, LF).

1	2	3	4	5
A1	A2	A3	CR	LF

### Response types

A1	A2	A3	Code			Meaning
A	0	0	41H	30H	30H	Successful completion
E	0	1	45H	30H	31H	* Command error (when an errant command is received)
E	0 - 9	0 - 9	45H	30H 5 39H	30H 5 39H	(Errors other than E01) * Numeric format error * Processing interrupted * Processing terminated abnormally * Other errors

### (2) ACK/NAK format

Consists of 1 character (Does not contain a terminator).

1
A1

### Response types

A1	Code	Meaning
ACK	06H	Successful completion
NAK	15H	* Command error (when an errant command is received) * Numeric format error * Processing interrupted * Processing terminated abnormally * Other errors

\* The following descriptions and tables assume that the A00/Exx format is selected as the response format.

Think of A00 as ACK, and E00 as NAK.



### 12.5.3 Command format

#### (1) Tare range (zero-setting) command

C1	C2	Code		Description	Value	Response
T	(SP)	54H	20H	<ul style="list-style-type: none"> <li>Set Tare Range</li> <li>Adjust the zero-point</li> </ul>	None	A00: Successful completion E01: Command error E04: A tare range cannot be set or the zero-point cannot be adjusted (because of a range violation or a weight error).

#### (2) Set output control command

C1	C2	Code		Description	
0	0	4FH	30H	Stop output	A00: Successful completion E01: Command error E02: Interval time error (only with OA or OB)
0	1	4FH	31H	Output continuously at all times	
0	2	4FH	32H	Output continuously if stable (Stop output if unstable)	
0	3	4FH	33H	Output once when the <u>Print</u> key is pressed (whether the balance is stable or unstable).	
0	4	4FH	34H	Output once when the balance is stable. After any sample is lifted causing the display to show zero or less value, outputs the next time any sample is put on again and the balance is stabilized.	
0	5	4FH	35H	Output once when the balance is stable. Stop output when unstable. Output once when the balance is stabilized again (the output includes zero) even if it is not reloaded.	
0	6	4FH	36H	Output once when the balance is stable. Output continuously when unstable. Output is stopped after a single output when the balance is stable even if it is not reloaded.	
0	7	4FH	37H	Output once when the <u>Print</u> key is pressed if the balance is stable.	
0	8	4FH	38H	Output once immediately	
0	9	4FH	39H	Output once after stabilized	
0	A	4FH	41H	Interval function (Output once each time the output time elapses.)	
0	B	4FH	42H	Interval function (Output once if stable each time the output time elapses.)	

- \* Commands O0 to O7 control output in the same way as the settings in Function 1.
- \* Commands O8 and O9 are used to request data from the balance.
- \* Once executed, O0 to O7, OA, or OB commands is held until the next command is input. However, output control is reset to the function setting when the balance is turned on again.
- \* When command OA or OB is input, the balance initiates the interval function, and when the same command is input again, it terminates the function.
- \* Once either an O8 or O9 command is executed, it returns to the state of "O0".

## 12.5.4 Set measurement mode command

Command Main Body				Description	Value	Response
1st character	2nd character	Code				
M	1	4DH	31H	Set to Mode 1	None	A00: Successful completion E01: Command error E02: (error)
M	2	4DH	32H	Set to Mode 2		
M	3	4DH	33H	Set to Mode 3		
M	4	4DH	34H	Set to Mode 4		

\* The measurement mode to be activated by the above mode settings 1 to 4 depends on the weighing mode currently in use.

Mode setting	Weighing machine	Parts counting	Percentage weighing	Unit converting	Gravimeter
Mode 1	Weight measuring	Weight measuring	Weight measuring	Weight measuring	(Error)
Mode 2	Gross weight display	Parts counting	Percent measuring	Coefficient multiplying	(Error)
Mode 3	Cumulative weight <sup>*1</sup>	Cumulative count display <sup>*1</sup>	Cumulative percent display <sup>*1</sup>	Cumulative sum <sup>*1</sup>	(Error)
Mode 4	Display in unit B <sup>*2</sup>	Average unit weight	(Error)	(Error)	(Error)

\*1 Mode 3 (M3) can be specified only when the addition function is to be used. If the addition function is not enabled, it results in an error.

\*2 If a unit is not specified for unit B, the balance is set for the weight measuring mode. When you specify a mode that is not supported by the current Weighing Mode, an error is returned from the balance.

## 12.5.5 Request date or time output command

Command Main Body				Description	Value	Response
1st character	2nd character	Code				
D	D	44H	44H	Date output request	None	A00: Successful completion E01: Command error
D	T	44H	54H	Time output request		

### (1) Date format

DATE:yyyy.mm.dd(CR)(LF)

English

㊦㊧㊨:yyyy.mm.dd(CR)(LF)

Japanese (Katakana)

### (2) Time data

TIME:(SP)(SP)(SP)(SP)(SP)hh:mm(CR)(LF)

English

㊦ ㊧㊨:(SP)(SP)(SP)(SP)(SP)hh:mm(CR)(LF)

Japanese (Katakana)

The above data is accompanied by a control command of our printer when output.

\*1 You can select whether you want to output data in English or in Japanese (katakana) in the Printed Language [*E* *J* *P.F.*] Function setting.

\*2 Date format depends on the setting of Date Display in Function 1 [*F* *d* *R* *t* *E*].

### 12.5.6 Setting intervals

Command Main Body				Description	Value	Response
1st character	2nd character	Code				
I	A	49H	41H	Interval time setup	Interval time	A00: Successful completion E01: Command error E02: Interval time error

- \* Upon setting the time interval, mark off between the entries of command, hour, minute, and second using a comma, "," between each of the entries; e.g., IA, hh, mm, ss.

### 12.5.7 Span adjustment/test command

Command Main Body				Description	Value	Response
1st character	2nd character	Code				
C	0	43H	30H	Disables <b>Cal</b> key and command inputs.	None	A00: Successful completion E01: Command error E02: Operation is disabled. E03: Cancelled by operation E04: Abnormal completion
C	1	43H	31H	Span adjustment with built-in weight (Auto-span adjustment)		
C	2	43H	32H	Span adjustment with built-in weight		
C	3	43H	33H	Span Adjustment to use external weights		
C	4	43H	34H	Span Test to use external weights		

- \* This command takes time because the balance sends back a response after an appropriate operation is completed.
- \* If the **Cal** key is disabled in Span Adjustment/Span Test of Function 1 [ $\overline{F}$   $\overline{R}$ ], a span adjustment and test command does not work.
- \* If the **Cal** key is disabled, however, the setting is reset to the function setting when the balance is turned on again or [ $\overline{F}$   $\overline{R}$ ] is changed directly. Span adjustment/test commands are also disabled.

### 12.5.8 Sample input commands

T(SP)(CR)(LF)	Set tare range or adjust to zero-point.
O1(CR)(LF)	Set the balance to continuous output.
O8(CR)(LF)	Output data (once immediately).
IA,01,30,00(CR)(LF)	Set an interval of 01 hrs: 30 min: 00 sec.
OA(CR)(LF)	Start the interval function.
DD(CR)(LF)	Output the date.
DT(CR)(LF)	Output the time.

# 13 Use Printers

## 13.1 Printer Setting

- (1) Use CSP-160 or CSP-240 with the balance.
- (2) Set proper print functions (print control) with the balance referring to the instruction manual for your printer.  
The factory default of our printer is manual printing (printer control).
- (3) Make the baud rate and other settings compatible between the balance and the printer.

## 13.2 Saving CAL and span test results

<p>1 Set Function 1.</p> <div style="display: flex; justify-content: space-around;"> <div style="border: 1px solid black; padding: 5px; width: 100px; text-align: center;">E GLP I</div> <div style="border: 1px solid black; padding: 5px; width: 100px; text-align: center;">E l o u t I</div> </div>	<p>Set the ISO/GLP/GMP compatible item to [E GLP I] and set the GLP compatible item for measurement data to [E l o u t I].</p> <p>Then, press the <b>Set</b> key to return to the weight display.</p>
<p>2 Performs span adjustment/span test.</p> <div style="border: 1px solid black; padding: 5px; width: 150px; text-align: center;">b u s y </div>	<p>During output, it indicates [b u s y ] or an error display and it looks like the balance is frozen. However, this is a normal operation. Wait until the output operation is completed.</p>

## 13.3 Measurement Data Output in Compliance with ISO/GLP/GMP

<p>1 Set Function 1.</p> <div style="display: flex; justify-content: space-around;"> <div style="border: 1px solid black; padding: 5px; width: 100px; text-align: center;">E GLP I</div> <div style="border: 1px solid black; padding: 5px; width: 100px; text-align: center;">E 2 o d I</div> </div>	<p>Set the ISO/GLP/GMP compatible item to [E GLP I] and set the GLP compatible item for measurement data to [E 2 o d I].</p> <p>Then, press the <b>Set</b> key to return to the weight display.</p>
<p>2 Header output</p> <div style="display: flex; align-items: center;"> <div style="border: 1px solid black; padding: 5px; width: 150px; text-align: center;">H E R d </div> </div> <p>Continuous press</p>	<p>Press and hold the <b>Print</b> key for a while.</p> <p>The display indicates [H E R d] and a header is output.</p>
<p>3 Measurement data output</p>	<p>Measurement data can be output at any time during measurements.</p>
<p>4 Footer output</p> <div style="display: flex; align-items: center;"> <div style="border: 1px solid black; padding: 5px; width: 150px; text-align: center;">F o o t </div> </div> <p>Continuous press</p>	<p>Press and hold the <b>Print</b> key for a while after the measurement is completed.</p> <p>The display indicates [F o o t] and a footer is output.</p>

### Caution

- (1) Upon printing calibration/span test results or ISO/GLP/GMP compatible outputs, date and time data is also printed.  
Check the date and time set in the balance before printing. (Refer to Section 10: Date and Time Setup in page 40).
- (2) Do not press any keys on the printer during printing.

\* Refer to Section 14 "Output in Compliance with ISO/GLP/GMP" described from page 59 for sample prints.

# 14 Output in Compliance with ISO/GLP/GMP

When a span adjustment or span test is successfully completed, the balance presents the results according to the following format. If span adjustment or a span test is not completed successfully, no data will be printed. The "X" entries in the following tables will vary with the model or the weight used. The line of external weight "ERR" or weight error is printed only when weight error is used according to the Function 2 setting.

## (1) Span adjustment with the built-in weight (auto-span adjustment)

English

1	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	
2	*	*	C	A	L	I	B	R	A	T	I	O	N	*	*	
3																
4	D	A	T	E	:	2	0	0	5	.	0	7	.	1	1	
5	T	I	M	E	:						1	6	:	5	6	
6	S	H	I	N	K	O					D	E	N	S	H	I
7	T	Y	P	E	:											
8	X	X	X	X	X	X	X	H	J	-	6	2	0	0	E	
9	S	/	N	:		X	X	X	X	X	X	X	X	X	X	
10	I	D	:							X	X	X	X	X	X	
11																
12	C	A	L	.	I	N	T	E	R	N	A	L				
13	R	E	F	:												
14			X	X	X	X	X	X	X	X	X	X	X	X	q	
15																
16	C	O	M	P	L	E	T	E								
17	D	A	T	E	:	2	0	0	5	.	0	7	.	1	1	
18	T	I	M	E	:						1	6	:	5	7	
19																
20	S	I	G	N	A	T	U	R	E							
21																
22																
23																
24	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	
25																
26																
27																
28																

Japanese (Katakana)

1	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	
2	*	*	*					コ	ウ	セ	イ			*	*	*
3																
4	ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1	
5	シ	"	コ	ク	:						1	6	:	5	6	
6			S	H	I	N	K	O			D	E	N	S	H	I
7	カ	タ	シ	キ	:											
8	X	X	X	X	X	X	X	H	J	-	6	2	0	0	E	
9	セ	イ	ハ	"	ン	X	X	X	X	X	X	X	X	X	X	
10	I	D	:								X	X	X	X	X	
11																
12	コ	ウ	セ	イ	(	カ	"	イ	フ	"	フ	ン	ト	"	ウ	)
13	キ	シ	"	ユ	ン	:										
14			X	X	X	X	X	X	X	X	X	X	X	X	X	q
15																
16	シ	ユ	ウ	リ	ヨ	ウ										
17	ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1	
18	シ	"	コ	ク	:						1	6	:	5	7	
19																
20	シ	ヨ	メ	イ												
21																
22																
23																
24	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*
25																
26																
27																
28																

## (2) Span adjustment with an external weight

English

1	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	
2	*	*	C	A	L	I	B	R	A	T	I	O	N	*	*	
3																
4	D	A	T	E	:	2	0	0	5	.	0	7	.	1	1	
5	T	I	M	E	:						1	6	:	5	6	
6	S	H	I	N	K	O					D	E	N	S	H	I
7	T	Y	P	E	:											
8	X	X	X	X	X	X	X	H	J	-	6	2	0	0	E	
9	S	/	N	:		X	X	X	X	X	X	X	X	X	X	
10	I	D	:							X	X	X	X	X	X	
11																
12	C	A	L	.	E	X	T	E	R	N	A	L				
13	R	E	F	:												
14			X	X	X	X	X	X	X	X	X	X	X	X	q	
15	E	R	R	:	X	X	X	X	X	3	.	2	1	m	q	
16																
17	C	O	M	P	L	E	T	E								
18	D	A	T	E	:	2	0	0	5	.	0	7	.	1	1	
19	T	I	M	E	:						1	6	:	5	6	
20																
21	S	I	G	N	A	T	U	R	E							
22																
23																
24																
25	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*
26																
27																
28																
29																

Japanese (Katakana)

1	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	
2	*	*	*					コ	ウ	セ	イ			*	*	*
3																
4	ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1	
5	シ	"	コ	ク	:						1	6	:	5	6	
6			S	H	I	N	K	O			D	E	N	S	H	I
7	カ	タ	シ	キ	:											
8	X	X	X	X	X	X	X	H	J	-	6	2	0	0	E	
9	セ	イ	ハ	"	ン	X	X	X	X	X	X	X	X	X	X	
10	I	D	:								X	X	X	X	X	
11																
12	コ	ウ	セ	イ	(	カ	"	イ	フ	"	フ	ン	ト	"	ウ	)
13	キ	シ	"	ユ	ン	:										
14			X	X	X	X	X	X	X	X	X	X	X	X	X	q
15	キ	サ	:	X	X	X	X	X	X	3	.	2	1	m	q	
16																
17	シ	ユ	ウ	リ	ヨ	ウ										
18	ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1	
19	シ	"	コ	ク	:						1	6	:	5	7	
20																
21	シ	ヨ	メ	イ												
22																
23																
24																
25	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*
26																
27																
28																
29																

(3) Span test with the built-in weight

English

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15
1															
2	*	*	C	A	L	.	T	E	S	T	*	*	*	*	*
3															
4	D	A	T	E	:	2	0	0	5	.	0	7	.	1	1
5	T	I	M	E	:						1	6	:	5	6
6	S	H	I	N	K	O	.	D	E	N	S	H	I		
7	T	Y	P	E	:										
8	X	X	X	X	X	X	X	H	J	-	6	2	0	0	E
9	S	/	N	:		X	X	X	X	X	X	X	X	X	X
10	I	D	:						X	X	X	X	X	X	X
11															
12	C	A	L	.	I	N	T	.	T	E	S	T			
13	R	E	F	:											
14			X	X	X	X	X	X	X	X	X	X	X		q
15	D	I	F	F	:										
16			X	X	X	X	X	X	X	X	X	X	X		q
17															
18	C	O	M	P	L	E	T	E							
19	D	A	T	E	:	2	0	0	5	.	0	7	.	1	1
20	T	I	M	E	:						1	6	:	5	7
21															
22	S	I	G	N	A	T	U	R	E						
23															
24															
25															
26	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*
27															
28															
29															
30															

Japanese (Katakana)

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	
1																
2	*	*	*						テ	ス	ト			*	*	*
3																
4	ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1	
5	シ	"	コ	ク	:							1	6	:	5	6
6	S	H	I	N	K	O	.	D	E	N	S	H	I			
7	カ	タ	シ	キ	:											
8	X	X	X	X	X	X	X	H	J	-	6	2	0	0	E	
9	セ	イ	ハ	"	ン	X	X	X	X	X	X	X	X	X	X	X
10	I	D	:							X	X	X	X	X	X	X
11																
12	テ	ス	ト	(	ナイ	フ	"	フ	ン	ト	"	ウ	)			
13	キ	シ	"	ユ	ン	:										
14			X	X	X	X	X	X	X	X	X	X	X	X		q
15	コ	"	サ	:												
16			X	X	X	X	X	X	X	X	X	X	X	X		q
17																
18	シ	ユ	ウ	リ	ョ	ウ										
19	ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1	
20	シ	"	コ	ク	:							1	6	:	5	7
21																
22	シ	ョ	メ	イ												
23																
24																
25																
26	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*
27																
28																
29																
30																

(4) Span test with an external weight

English

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	
1																
2	*	*	C	A	L	.	T	E	S	T	*	*	*	*	*	
3																
4	D	A	T	E	:	2	0	0	5	.	0	7	.	1	1	
5	T	I	M	E	:						1	6	:	5	6	
6	S	H	I	N	K	O	.	D	E	N	S	H	I			
7	T	Y	P	E	:											
8	X	X	X	X	X	X	X	H	J	-	6	2	0	0	E	
9	S	/	N	:		X	X	X	X	X	X	X	X	X	X	
10	I	D	:						X	X	X	X	X	X	X	
11																
12	C	A	L	.	E	X	T	.	T	E	S	T				
13	R	E	F	:												
14			X	X	X	X	X	X	X	X	X	X	X		q	
15	E	R	R	:	X	X	X	X	X	3	.	2	1	m	q	
16	D	I	F	F	:											
17			X	X	X	X	X	X	X	X	X	X	X		q	
18																
19	C	O	M	P	L	E	T	E								
20	D	A	T	E	:	2	0	0	5	.	0	7	.	1	1	
21	T	I	M	E	:						1	6	:	5	7	
22																
23	S	I	G	N	A	T	U	R	E							
24																
25																
26																
27	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*
28																
29																
30																
31																

Japanese (Katakana)

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	
1																
2	*	*	*						テ	ス	ト			*	*	*
3																
4	ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1	
5	シ	"	コ	ク	:							1	6	:	5	6
6	S	H	I	N	K	O	.	D	E	N	S	H	I			
7	カ	タ	シ	キ	:											
8	X	X	X	X	X	X	X	H	J	-	6	2	0	0	E	
9	セ	イ	ハ	"	ン	X	X	X	X	X	X	X	X	X	X	X
10	I	D	:							X	X	X	X	X	X	X
11																
12	テ	ス	ト	(	カ	"	イ	フ	"	フ	ン	ト	"	ウ	)	
13	キ	シ	"	ユ	ン	:										
14			X	X	X	X	X	X	X	X	X	X	X	X		q
15	キ	サ	:	X	X	X	X	X	3	.	2	1	m	q		
16	コ	"	サ	:												
17			X	X	X	X	X	X	X	X	X	X	X		q	
18																
19	シ	ユ	ウ	リ	ョ	ウ										
20	ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1	
21	シ	"	コ	ク	:							1	6	:	5	7
22																
23	シ	ョ	メ	イ												
24																
25																
26																
27	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*
28																
29																
30																
31																

(5) Printing of measurement data

1) Header

English

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15		
1																	
2			S	H	I	N	K	O		D	E	N	S	H	I		
3			T	Y	P	E	:	X	X	H	J	-	6	2	0	0	E
4			S	/	N	:		X	X	X	X	X	X	X	X	X	X
5			I	D	:					X	X	X	X	X	X	X	X
6																	
7			S	T	A	R	T										
8			D	A	T	E	:	2	0	0	5	.	0	7	.	1	1
9			T	I	M	E	:						1	6	:	5	5
10																	

Japanese (Katakana)

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15		
1																	
2			S	H	I	N	K	O		D	E	N	S	H	I		
3			カ	タ	シ	キ	:	X	X	H	J	-	6	2	0	0	E
4			セ	イ	ハ	"	ン	X	X	X	X	X	X	X	X	X	X
5			イ	ド	:					X	X	X	X	X	X	X	X
6																	
7			カ	イ	シ												
8			ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1
9			シ	"	コ	ク	:						1	6	:	5	5
10																	

2) Footer

English

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15		
1																	
2			E	N	D												
3			D	A	T	E	:	2	0	0	5	.	0	7	.	1	1
4			T	I	M	E	:					1	6	:	5	6	
5																	
6			S	I	G	N	A	T	U	R	E						
7																	
8																	
9																	
10			*	*	*	*	*	*	*	*	*	*	*	*	*	*	*
11																	
12																	
13																	
14																	

Japanese (Katakana)

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15		
1																	
2			シ	ユ	ウ	リ	ョ	ウ									
3			ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1
4			シ	"	コ	ク	:					1	6	:	5	6	
5																	
6			シ	ョ	メ	イ											
7																	
8																	
9																	
10			*	*	*	*	*	*	*	*	*	*	*	*	*	*	*
11																	
12																	
13																	
14																	

(6) Calibration of Built-in Weight

English

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15		
1																	
2			*	*	*	*	R	E	F	.	C	A	L	*	*	*	*
3																	
4			D	A	T	E	:	2	0	0	5	.	0	7	.	1	1
5			T	I	M	E	:					1	6	:	5	6	
6			S	H	I	N	K	O		D	E	N	S	H	I		
7			T	Y	P	E	:	X	X	H	J	-	6	2	0	0	E
8			X	X	X	X	X	X	X	X	X	X	X	X	X	X	X
9			S	/	N	:		X	X	X	X	X	X	X	X	X	X
10			I	D	:					X	X	X	X	X	X	X	X
11																	
12			R	E	F	:											
13						X	X	X	X	X	X	X	X	X			g
14			E	R	R	:	X	X	X	X	X	3	.	2	1	m	g
15																	
16			C	O	M	P	L	E	T	E							
17			D	A	T	E	:	2	0	0	5	.	0	7	.	1	1
18			T	I	M	E	:					1	6	:	5	8	
19																	
20			S	I	G	N	A	T	U	R	E						
21																	
22																	
23																	
24			*	*	*	*	*	*	*	*	*	*	*	*	*	*	*
25																	
26																	
27																	
28																	

Japanese (Katakana)

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15		
1																	
2			*	ナ	イ	フ	"	フ	ン	ト	"	ウ	コ	ウ	セ	イ	*
3																	
4			ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1
5			シ	"	コ	ク	:					1	6	:	5	6	
6			S	H	I	N	K	O		D	E	N	S	H	I		
7			カ	タ	シ	キ	:	X	X	H	J	-	6	2	0	0	E
8			X	X	X	X	X	X	X	X	X	X	X	X	X	X	X
9			セ	イ	ハ	"	ン	X	X	X	X	X	X	X	X	X	X
10			イ	ド	:					X	X	X	X	X	X	X	X
11																	
12			キ	シ	"	ユ	ン	:									
13					X	X	X	X	X	X	X	X	X	X			g
14			キ	サ	:	X	X	X	X	X	X	3	.	2	1	m	g
15																	
16			シ	ユ	ウ	リ	ョ	ウ									
17			ヒ	ツ	"	ケ	:	2	0	0	5	.	0	7	.	1	1
18			シ	"	コ	ク	:					1	6	:	5	8	
19																	
20			シ	ョ	メ	イ											
21																	
22																	
23																	
24			*	*	*	*	*	*	*	*	*	*	*	*	*	*	*
25																	
26																	
27																	
28																	

## 15 Operate on Batteries



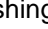

---

This function can only be used with a balance equipped with optional batteries. The battery option is unavailable for HJR series units.

### 15.1 Specifications

- Built-in batteries
- Charge time: About 15 hours
- Operation time: About 6 hours of continuous operation
- Can be recharged: More than 300 times

### 15.2 Recharging Batteries

- A [  ] icon is displayed when the balance is operated on batteries. As the battery gets weaker, the sign is changed into [  ] or [  ] (flashing). When the icon changes to [  ] (flashing), follow the steps below to charge the batteries.
  - (1) Connect the AC adapter supplied with the balance.
  - (2) Turn the balance off.
  - (3) Wait about 15 hours until the batteries are fully recharged.

### 15.3 Precautions Relating to Use

1. Disconnect the AC adapter after the batteries are recharged. The balance operated on batteries is slightly charged while it is in use (turned on). Do not continue to charge batteries after they are recharged, because that may shorten their life.
2. When you use the balance for the first time after purchase, battery operation time may be shorter than usual because of natural discharge.
3. If nothing is displayed or the display is turned off in a minute after you turn on the balance or the display is turned off after a succession of short beeps, the batteries may be weak. Recharge the batteries or operate the balance on the AC adapter.



#### **Caution** Follow the following instructions for safe use:

1. Do not disassemble or modify the batteries. Avoid the reverse wiring of positive and negative poles or a short circuit. Such mishandling may damage the batteries, or cause the balance to fail or ignite.
2. Be sure to use the AC adapter supplied with the balance. A different AC adapter may cause the batteries to generate heat or explode.
3. Do not put the batteries into a fire because they may explode.



## 16 Troubleshooting

\* Parentheses contain a page to refer to.

Symptom	Cause	Measures to Take
No display	<ul style="list-style-type: none"> <li>The AC adapter is not connected.</li> <li>The batteries are exhausted.</li> </ul>	<p>→ Ensure that the AC adapter is connected.</p> <p>→ Recharge the batteries (page 62).</p>
<p>The display is unstable.</p> <p>The display takes time to be stabilized.</p> <p>The display does not move with the M sign flashing.</p>	<ul style="list-style-type: none"> <li>The balance is influenced by wind or vibration.</li> <li>The balance is placed on an unstable base.</li> <li>The pan, tare, or sample put on the balance is in contact with an external object.</li> </ul>	<p>→ Read Precautions Relating to Use (pages 2 to 4).</p>
Weight indication contains an error.	<ul style="list-style-type: none"> <li>Tare range is set incorrectly.</li> <li>The adjusters are not settled, and the balance is not kept horizontal.</li> <li>The display error is caused because the balance has not been used for a long period of time or has been relocated to another location.</li> </ul>	<p>→ Set the tare range correctly.</p> <p>→ Check that the balance is kept horizontal (page 8).</p> <p>→ Perform span adjustment of the balance (page 37).</p>
[ $\square$ - E r r ] appears before the capacity is reached.	<ul style="list-style-type: none"> <li>The capacity is reached by the weight including the tare weight. (Measurement range is for the sum of container and sample weight.)</li> <li>The mechanism section was damaged for some reason.</li> </ul>	<p>→ Check the total weight.</p> <p>→ Review the container.</p>
[ $\square$ - E r r ] is displayed.	<ul style="list-style-type: none"> <li>There is a foreign object trapped between the pan (or pan base) and the balance.</li> <li>The mechanism section was damaged for some reason.</li> </ul>	<p>→ Remove the pan and examine the surface on which the pan was placed.</p>
[ $\square$ - E r r ] is displayed.	<ul style="list-style-type: none"> <li>The internal clock has failed.</li> </ul>	<p>→ Contact our technical service personnel or your local dealer.</p>
<p>[<math>\square</math> - E r r ] is displayed.</p> <p>[<math>\square</math> - E r r ] is displayed.</p>	<ul style="list-style-type: none"> <li>The balance is influenced by static electricity or noise.</li> <li>The electric components of the balance have failed.</li> </ul>	<p>→ Contact our technical service personnel or your local dealer.</p>
<p>The display is turned off when the balance is operated on batteries.</p> <p>The [<math>\square</math>] icon blinks.</p> <p>No display</p>	<ul style="list-style-type: none"> <li>The display was turned off by the auto power off function.</li> <li>The batteries were exhausted.</li> </ul>	<p>→ Turn the power on again. Disable the auto power off function if desired (page 41).</p> <p>→ Recharge the batteries (page 62).</p> <p>→ Operate the balance on the AC adapter.</p>
[ $\square$ - E r r ] is displayed.	The built-in weight drive unit has have failed.	<p>→ Contact your local dealer.</p>

# 17 Specifications

## 17.1 Basic Specifications

Model \ Item	HJ-220E	HJ-320E	HJ-420E	HJ-620E	HJ-1200E	HJ-2200E	HJ-3200E	HJ-4200E
	HJR-220E	HJR-320E	HJR-420E	HJR-620E	HJR-1200E	HJR-2200E	HJR-3200E	HJR-4200E
Weighing capacity [g]	220g	320g	420g	620g	1200g	2200g	3200	4200
Minimum readability [g]	0.001g	0.001g	0.001g	0.001g	0.01g	0.01g	0.01g	0.01g
Parts counting Minimum unit weight [g]	0.001g	0.001g	0.001g	0.001g	0.01g	0.01g	0.01g	0.01g
Percentage Lower weight limit [g]	0.1g	0.1g	0.1g	0.1g	1g	1g	1g	1g
Size of measurement pan	120 × 140[mm]				200 × 200[mm]			
Span adjustment method	HJ series: External weight is used HJR series: Internal or external weight is used							

Model \ Item	HJ-6200E	HJ-8200E	HJ-10KE	HJ-12KE	HJ-15KE	HJ-21KE	HJ-31KE
	Weighing capacity [g]	6200g	8200g	10000g	12000g	15000g	21000g
Minimum readability [g]	0.01g	0.1g	0.05g	0.1g	0.1g	0.1g	0.1g
Parts counting Minimum unit weight [g]	0.01g	0.1g	0.05g	0.1g	0.1g	0.1g	0.1g
Percentage Lower weight limit [g]	1g	10g	5g	10g	10g	10g	10g
Size of measurement pan	200 × 200[mm]					220 × 250[mm]	
Span adjustment method	Span adjustment with external weight						

## 17.2 Common Specifications

- (1) Weight measurement method: ..... Tuning fork system
- (2) Tare range ..... Full capacity
- (3) Display: ..... The FIP segment can display up to 8 digits.  
The segment is 12.5 mm high, and displays 8 digits for a weight.  
Can display a bar graph of up to 40 bars and various messages.
- (4) Overload display: ..... [  $\square$  -  $\bar{E}$  r r ] is displayed when the weighing capacity is exceeded by 9 divisions.
- (5) Output: ..... Bi-direction RS232C output, output for peripherals
- (6) Supported printer: ..... CSP-160, CSP-240
- (7) Operation temperature/humidity range: -5 to 35°C, 80%RH or lower
- (8) AC adapter: ..... Dedicated AC adapter: 12VDC/120VAC or  
12VDC/230VAC
- (9) Options .....
  - Hook for hanging \*1
  - Battery option \*2
  - Buzzer option
  - Limit contact output option
  - Full pack option
  - RS422 option \*3

\*1 For 1200g ~ 15kg Capacity models only.

(Inside is installed hook for under weighing for 220g ~ 620g, 21kg and 31kg model.)

\*2 unavailable for HJR series

\*3 When the RS422 option is selected, the D-SUB9P connector is replaced with an RS422 connector.

## 17.3 Minimum Display by Unit of Measurement

Model Unit_of measurement	HJ- 220E	HJ- 320E	HJ- 420E	HJ- 620E	HJ- 1200E	HJ- 2200E	HJ- 3200E	HJ- 4200E
	HJR- 220E	HJR- 320E	HJR- 420E	HJR- 620E	HJR- 1200E	HJR- 2200E	HJR- 3200E	HJR- 4200E
<b>g</b>	220 0.001	320 0.001	420 0.001	620 0.001	1200 0.01	2200 0.01	3200 0.01	4200 0.01
<b>kg</b>	0.22 0.000001	0.32 0.000001	0.42 0.000001	0.62 0.000001	1.2 0.00001	2.2 0.00001	3.2 0.00001	4.2 0.00001
<b>ct</b> (ct)	1100 0.005	1600 0.005	2100 0.005	3100 0.005	6000 0.05	11000 0.05	16000 0.05	21000 0.05
<b>oz</b> (oz)	7.7 0.00005	11 0.00005	14 0.00005	21 0.00005	42 0.00005	77 0.00005	110 0.00005	140 0.00005
<b>lb</b> (lb)	0.48 0.000005	0.7 0.000005	0.92 0.000005	1.3 0.000005	2.6 0.00005	4.8 0.00005	7 0.00005	9.2 0.00005
<b>oz t</b> (ozt)	7 0.00005	10 0.00005	13 0.00005	19 0.00005	38 0.00005	70 0.00005	100 0.00005	130 0.00005
<b>dwt</b> (dwt)	140 0.001	200 0.001	270 0.001	390 0.001	770 0.01	1400 0.01	2000 0.01	2700 0.01
<b>GN</b> (grain)	3300 0.02	4900 0.02	6400 0.02	9500 0.02	18000 0.2	33000 0.2	49000 0.2	64000 0.2
<b>ti</b> (Hong Kong)	5.8 0.00005	8.5 0.00005	11 0.00005	16 0.00005	32 0.00005	58 0.00005	85 0.00005	110 0.00005
<b>ti</b> (Singapore, Malaysia)	5.8 0.00005	8.4 0.00005	11 0.00005	16 0.00005	31 0.00005	58 0.00005	84 0.00005	110 0.00005
<b>ti</b> (Taiwan)	5.8 0.00005	8.5 0.00005	11 0.00005	16 0.00005	32 0.00005	58 0.00005	85 0.00005	110 0.00005
mom (momme)	58 0.0005	85 0.0005	110 0.0005	160 0.0005	320 0.005	580 0.005	850 0.005	1100 0.005
<b>to</b> (to)	18 0.0001	27 0.0001	36 0.0001	53 0.0001	100 0.001	180 0.001	270 0.001	360 0.001

\* The view of the table

Upper cell: Capacity
Lower cell: Readability

Model Unit_of measurament	HJ- 6200E	HJ- 8200E	HJ- 10KE	HJ- 12KE	HJ- 15KE	HJ- 21KE	HJ- 31KE
	<b>g</b>	6200 0.01	8200 0.1	10000 0.05	12000 0.1	15000 0.1	21000 0.1
<b>kg</b>	6.2 0.00001	8.2 0.0001	10 0.00005	12 0.0001	15 0.0001	21 0.0001	31 0.0001
<b>ct</b> (ct)	31000 0.05	41000 0.5	50000 0.5	60000 0.5	75000 0.5	100000 0.5	150000 0.5
<b>oz</b> (oz)	210 0.0005	280 0.005	350 0.002	420 0.005	520 0.005	740 0.005	1000 0.005
<b>lb</b> (lb)	13 0.00005	18 0.0005	22 0.0002	26 0.0005	33 0.0005	46 0.0005	68 0.0005
<b>oz t</b> (ozt)	190 0.0005	260 0.005	320 0.002	380 0.005	480 0.005	670 0.005	990 0.005
<b>dwt</b> (dwt)	3900 0.01	5200 0.1	6400 0.05	7700 0.1	9600 0.1	13000 0.1	19000 0.1
<b>GN</b> (grain)	95000 0.2	120000 2	150000 1	180000 2	230000 2	320000 2	470000 2
<b>tl</b> (Hong Kong)	160 0.0005	210 0.005	260 0.002	320 0.005	400 0.005	560 0.005	820 0.005
<b>tl</b> (Singapore, Malaysia)	160 0.0005	210 0.005	260 0.002	310 0.005	390 0.005	550 0.005	820 0.005
<b>tl</b> (Taiwan)	160 0.0005	210 0.005	260 0.002	320 0.005	400 0.005	560 0.005	820 0.005
mom (momme)	1600 0.005	2100 0.05	2600 0.02	3200 0.05	4000 0.05	5600 0.05	8200 0.05
<b>to</b> (to)	530 0.001	700 0.01	850 0.005	1000 0.01	1200 0.01	1800 0.01	2600 0.01

\* The view of the table

Upper cell: Capacity
Lower cell: Readability

## 18 Conversion Table of Units

Unit	Gram	carat	Ounce	Pound	troy ounce	Penny Weight
1g	1	5	0.03527	0.00220	0.03215	0.64301
1ct	0.2	1	0.00705	0.00044	0.00643	0.12860
1oz	28.34952	141.74762	1	0.06250	0.91146	18.22917
1lb	453.59237	2267.96185	16	1	14.58333	291.66667
1ozt	31.10348	155.51738	1.09714	0.06857	1	20
1dwt	1.55517	7.77587	0.05486	0.00343	0.05	1
1GN	0.06480	0.32399	0.00229	0.00014	0.00208	0.04167
1tl (HK)	37.429	187.145	1.32027	0.08252	1.20337	24.06741
1tl (SGP,Mal)	37.79936	188.99682	1.33333	0.08333	1.21528	24.30556
1tl (Taiwan)	37.5	187.5	1.32277	0.08267	1.20565	24.11306
1mom	3.75	18.75	0.13228	0.00827	0.12057	2.41131
1to	11.66380	58.31902	0.41143	0.02571	0.37500	7.5

unit	Grain	tael (Hong Kong)	tael (Singapore, Malaysia)	tael (Taiwan)	momme	Tola
1g	15.43236	0.02672	0.02646	0.02667	0.26667	0.08574
1ct	3.08647	0.00534	0.00529	0.00533	0.05333	0.01715
1oz	437.5	0.75742	0.75	0.75599	7.55987	2.43056
1lb	7000	12.11874	12	12.09580	120.95797	38.88889
1ozt	480	0.83100	0.82286	0.82943	8.29426	2.66667
1dwt	24	0.04155	0.04114	0.04147	0.41471	0.13333
1GN	1	0.00173	0.00171	0.00173	0.01728	0.00556
1tl (HK)	577.61774	1	0.99020	0.99811	9.98107	3.20899
1tl (SGP,Mal)	583.33333	1.00990	1	1.00798	10.07983	3.24074
1tl (Taiwan)	578.71344	1.00190	0.99208	1	10	3.21507
1mom	57.87134	0.10019	0.09921	0.1	1	0.32151
1to	180	0.31162	0.30857	0.31103	3.11035	1